


DYNAMIC ABUTMENT SOLUTIONS

GUIDED DASURGICAL KIT

WORK OFFSETS BY COMPATIBILITY



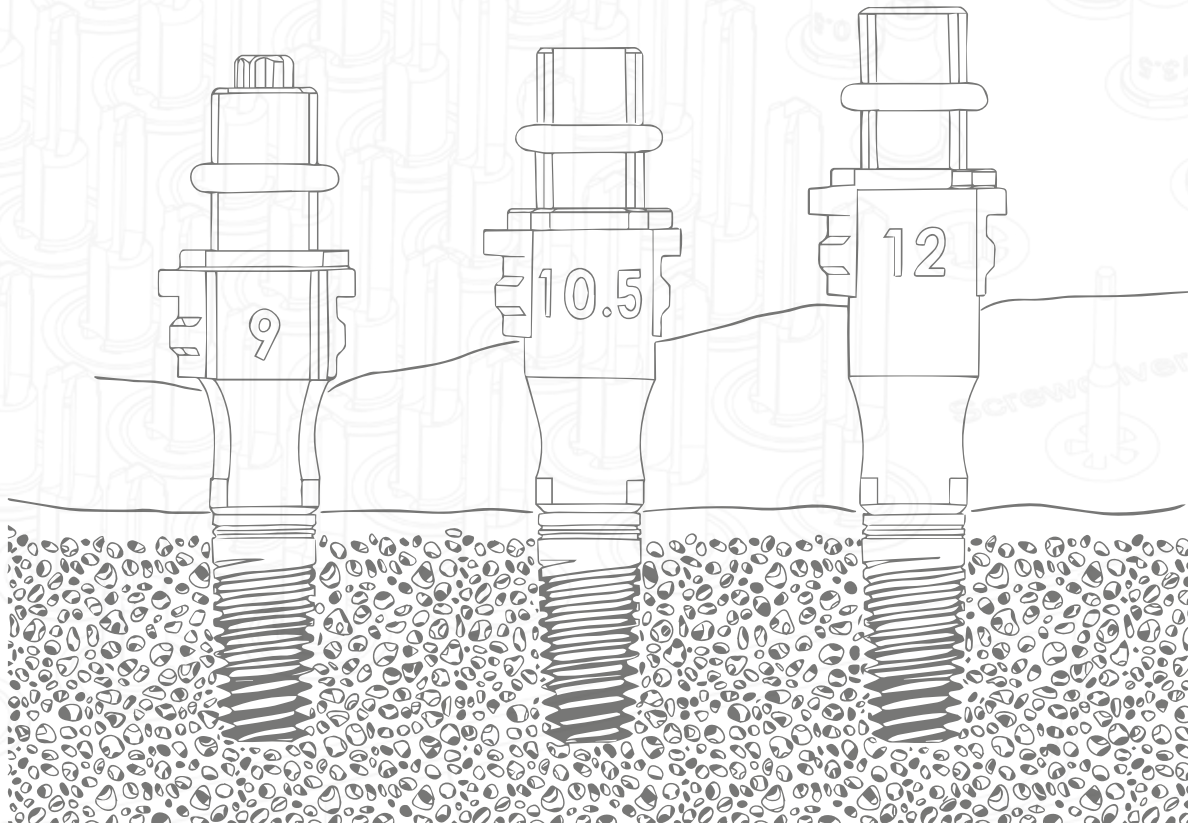

DYNAMIC ABUTMENT® SOLUTIONS

VIRGINIA WOOLF, 17 · 25005 · LLEIDA (SPAIN)

INTERNATIONAL +34 873 450 709

das@dynamicabutment.com

SPAIN +34 973 289 580

spain@dynamicabutment.com

Dynamic Abutment® Solutions is a trademark of the manufacturer Talladium España S.L

The information included in this catalogue is exclusively addressed to professionals in dental sector.

All commercial trademarks mentioned herein are fully registered by their respective companies, and the images that appear are just to provide an orientation.

Talladium reserves the right to modify this document without previous notice.

Talladium is not responsible for the inadequate execution of these products if the warning indications corresponding to every reference are not contemplated.

Request information on sales license and certificates according to the country of sale: das@dynamicabutment.com

INDEX

3	<i>ALPHABIO INTERNAL HEX CONNECTION (IH) - NEO/MULTINEO/SPIRAL/DFI/ATID/ICE</i>
11	<i>CORTEX INTERNAL HEX - DYNAMIX/CLASSIX</i>
17	<i>DIO - UF II REGULAR</i>
27	<i>HIOSEN - ET SYSTEM</i>
35	<i>MIS INTERNAL HEXAGON - SEVEN/M4</i>
41	<i>NEOBIOTECH - IS SYSTEM</i>
47	<i>NOBEL BIO CARE - ACTIVE/PARALLEL CONICAL CONNECTION/REPLACE CONICAL CONNECTION</i>
57	<i>OSSTEM - TS SYSTEM</i>
63	<i>PAL TOP IH - ADVANCED/ADVANCED +/DYNAMIC/ACTIVE</i>

REMARK: THERE ARE HERE SOME COMPATIBILITIES AS EXAMPLE.
Ask for other connections through e-mail das@dynamicabutment.com


ALPHABIO INTERNAL HEX CONNECTION (IH) - NEO/MULTINEO/SPIRAL/DFI/ATID/ICE

Implant Length (mm)	Mount (mm)	Depth (mm)	Drill Length (mm)
6	9	0	6
8	9.5	0 - 0.5	8.5
	12.5	0 - 3	11.5
10	9	0	10
	10.5	0 - 1.5	11.5
	12	0 - 3	13
11.5	9	0	11.5
	10.5	0 - 1.5	13
	13.5	0 - 3	16
13	9	0	13
	12	0 - 3	16
16	9	0	16
	12	0 - 3	19

Alphabio Internal Hex Connection (IH)		
Model	Ø Implant (mm)	Implant Length (mm)
NEO/MULTINEO	3.75	8/10/11.5/13/16
	4.2	8/10/11.5/13/16
SPIRAL	3.3	8/10/11.5/13/16
	3.75	8/10/11.5/13/16
	4.2	8/10/11.5/13/16
DFI	3.3	8/10/11.5/13/16
	3.75	8/10/11.5/13/16
	4.2	8/10/11.5/13/16
ATID	3.3	8/10/11.5/13/16
	3.75	8/10/11.5/13/16
	4.2	8/10/11.5/13/16
ICE	3.7N	10/11.5/13
	3.75	8/10/11.5/13/16
	4.2	6/8/10/11.5/13/16
	4.65	6/8/10/11.5/13

DAS_0040_B2


IMPLANT LENGTH 6 mm

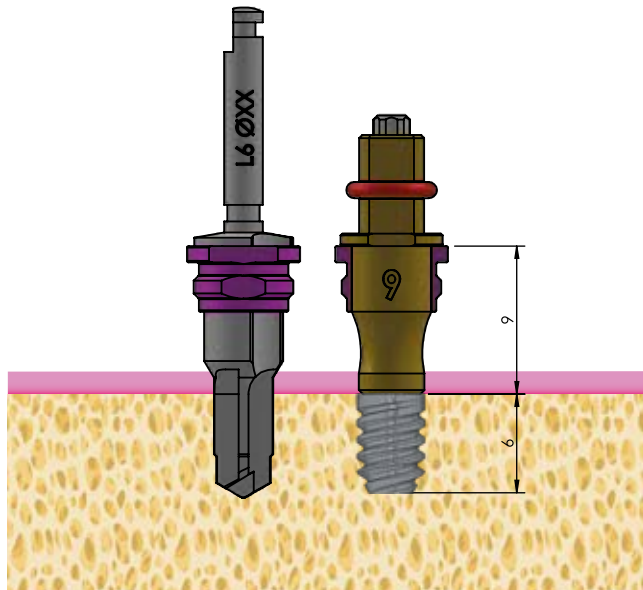


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
6	9	0	6

 $\varnothing XX = \varnothing \text{ Drill}$

IMPLANT LENGTH 8 mm

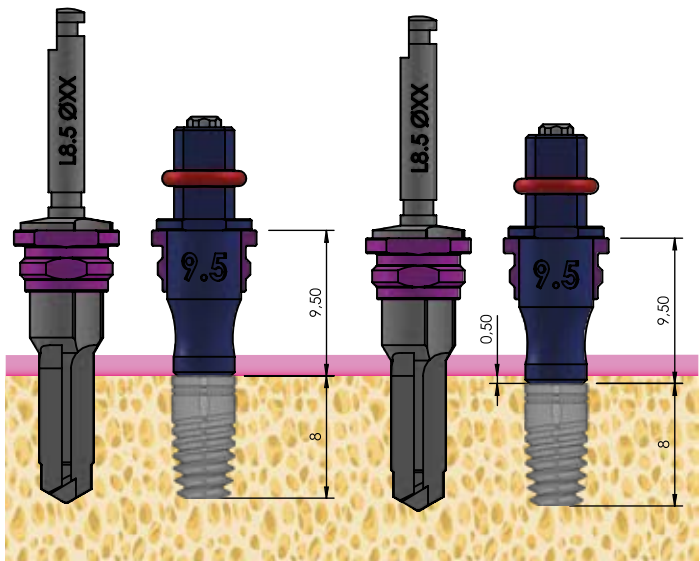


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8	9.5	0 - 0.5	8.5

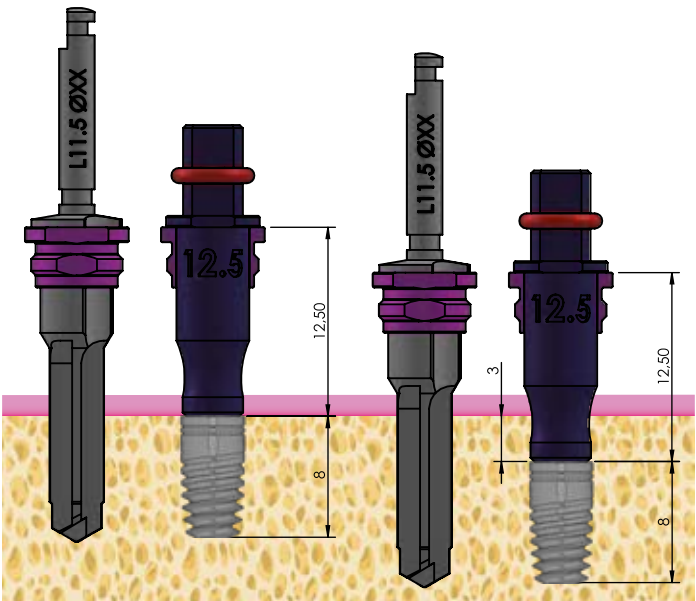


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8	12.5	0 - 3	11.5

Ø XX = Ø Drill


IMPLANT LENGTH 10 mm

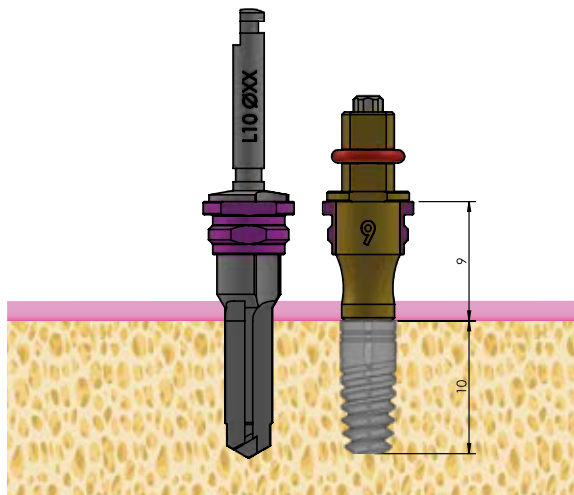

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	9	0	10

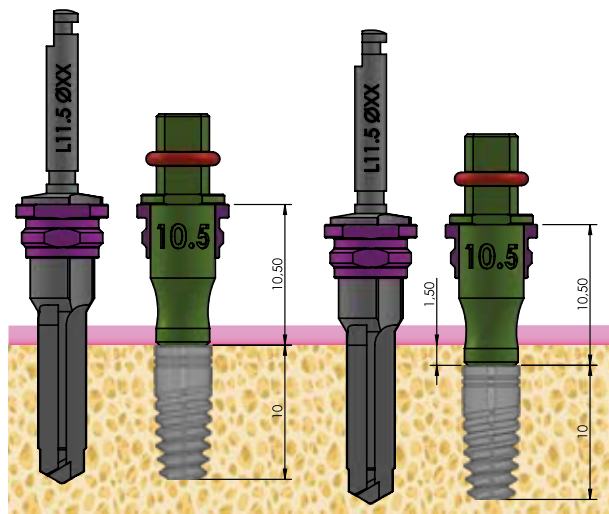


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	10.5	0 - 1.5	11.5

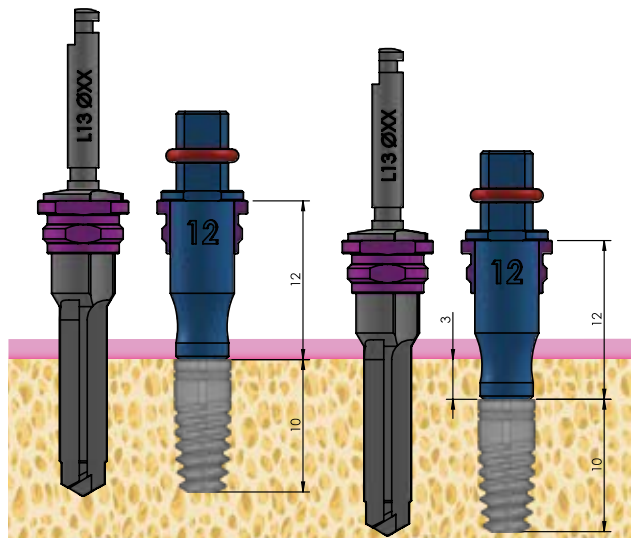


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	12	0 - 3	13

Ø XX = Ø Drill


IMPLANT LENGTH 11.5 mm

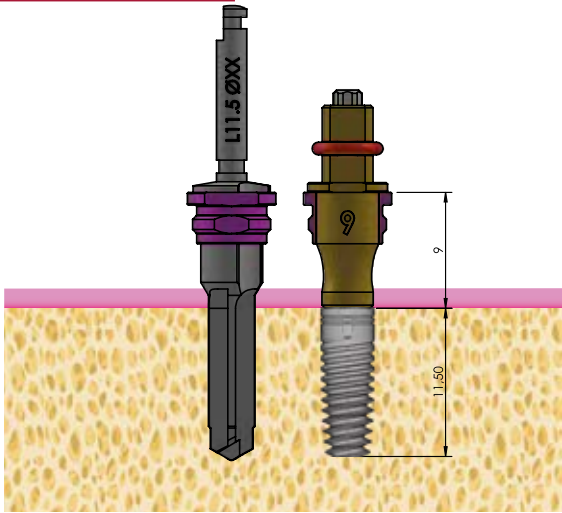

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	9	0	11.5

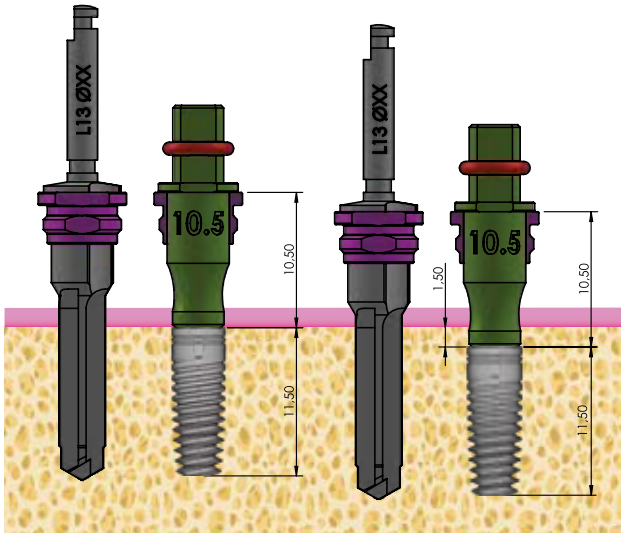


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	10.5	0 - 1.5	13

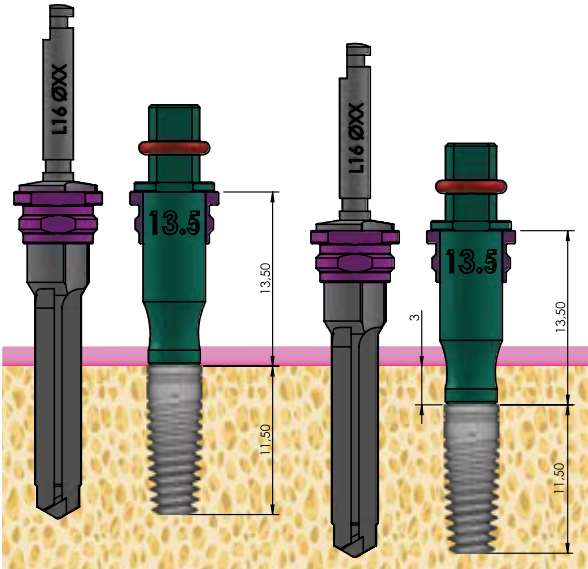


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	13.5	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 13 mm

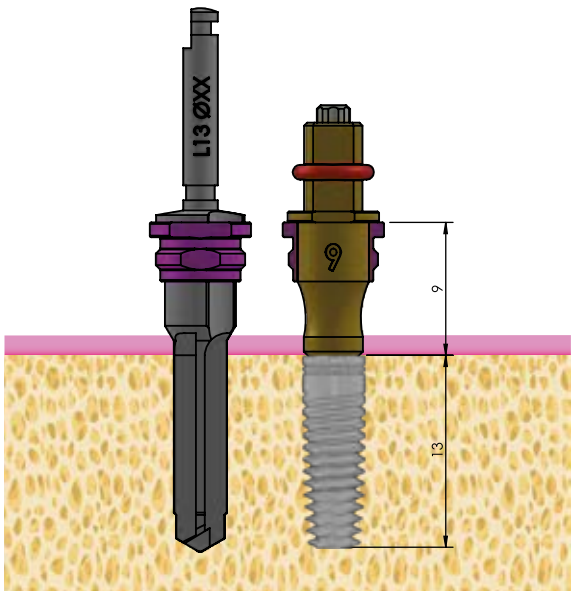


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	9	0	13

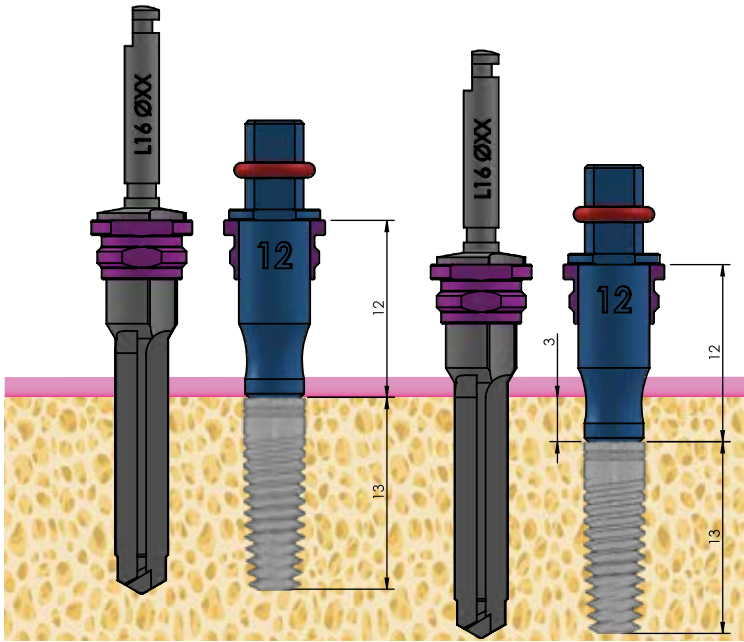


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	12	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 16 mm

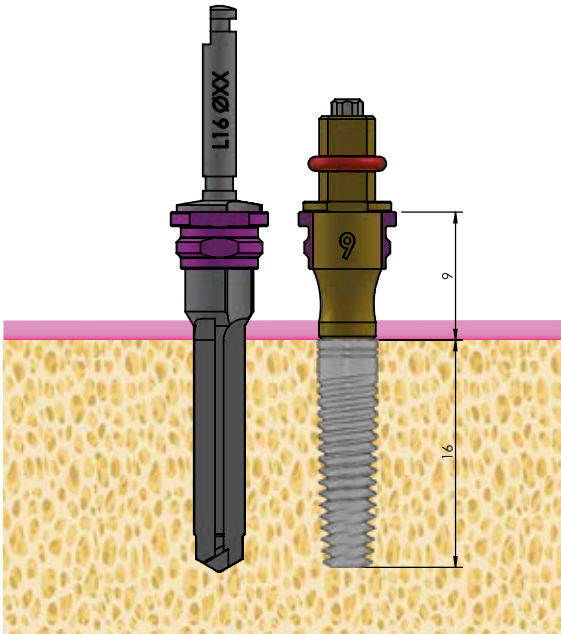

Implant Length (mm)	Mounter (mm)	Depth(mm)	Drill Length (mm)
16	9	0	16

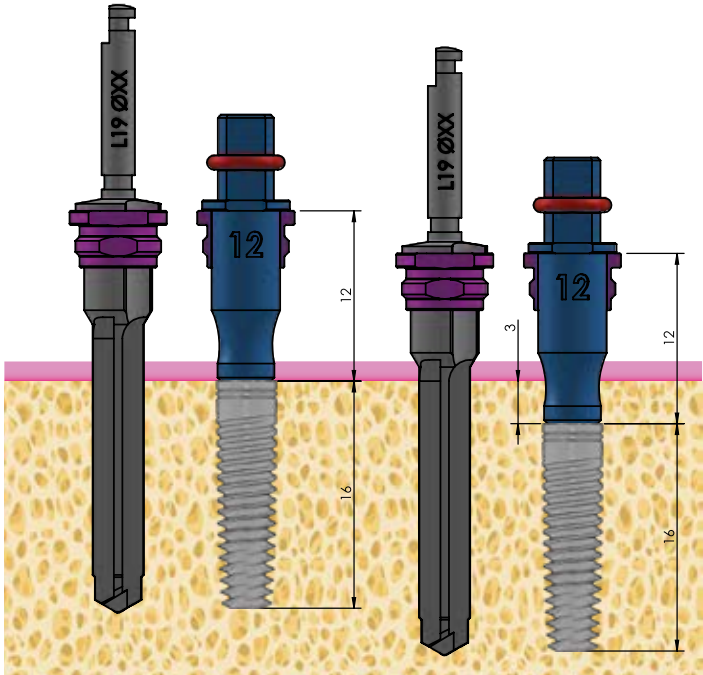

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
16	12	0 - 3	19

Ø XX = Ø Drill

NOTES


[illegible]

CORTEX INTERNAL HEX - DYNAMIX/CLASSIX

Implant Length (mm)	Mount (mm)	Depth (mm)	Drill Length (mm)
6	9	0	6
8	9.5	0 - 0.5	8.5
10	9	0	10
	10.5	0 - 1.5	11.5
	12	0 - 3	13
11.5	9	0	11.5
	10.5	0 - 1.5	13
	13.5	0 - 3	16
13	9	0	13
	12	0 - 3	16
16	9	0	16
	12	0 - 3	19


Cortex		
Model	Ø Implant (mm)	Implant Length (mm)
Dynamix	3.3	10/11.5/13/16
	3.8	8/10/11.5/13/16
	4.2	6/8/10/11.5/13/16
Classix	3.3	10/11.5/13/16
	3.8	8/10/11.5/13/16
	4.2	6/8/10/11.5/13/16

IMPLANT LENGTH 6 mm

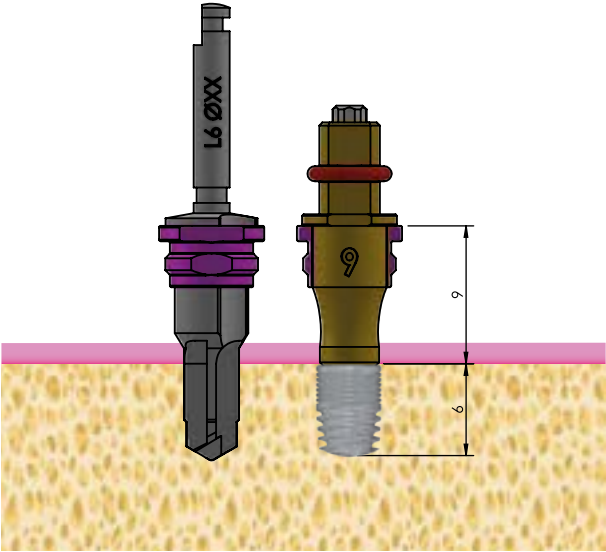

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
6	9	0	6


IMPLANT LENGTH 8 mm

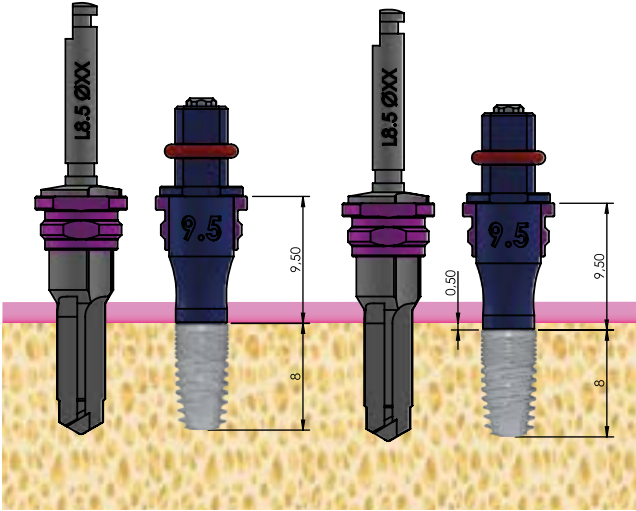


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8	9.5	0 - 0.5	8.5

Ø XX = Ø Drill


IMPLANT LENGTH 10 mm

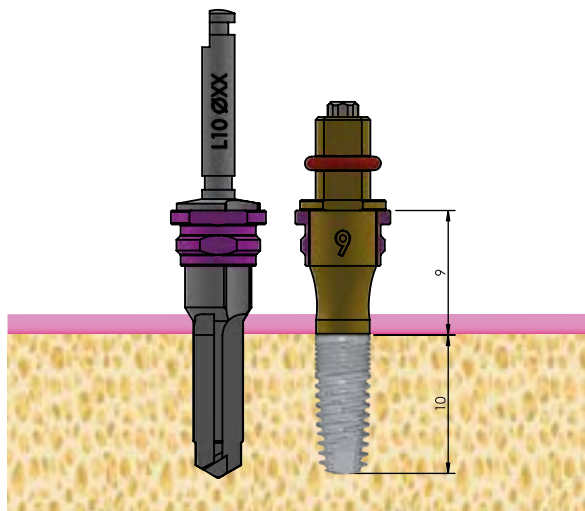

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	9	0	10

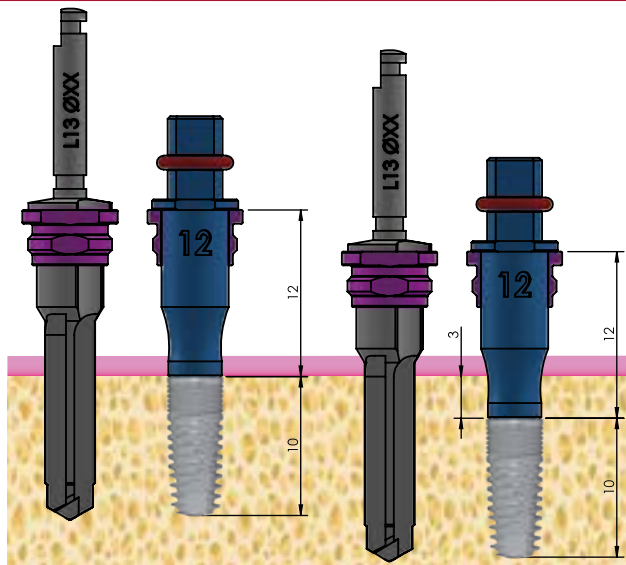


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	12	0 - 3	13

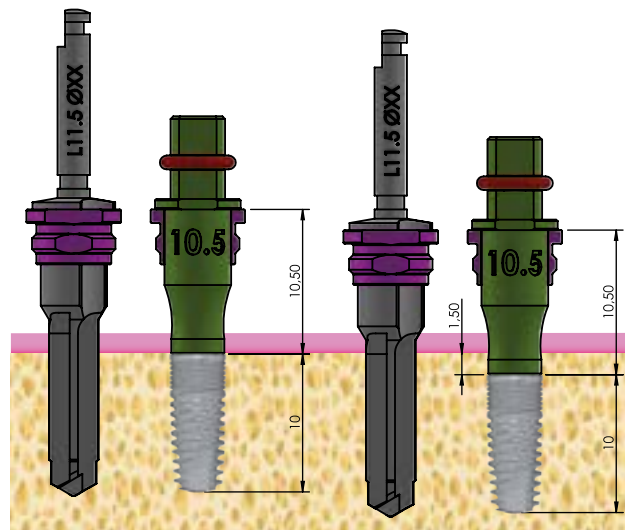


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	10.5	0 - 1.5	11.5

 $\emptyset XX = \emptyset \text{ Drill}$


IMPLANT LENGTH 11.5 mm

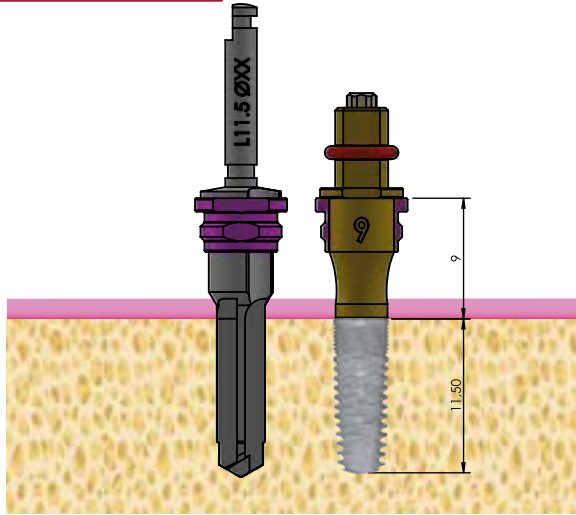

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	9	0	11.5

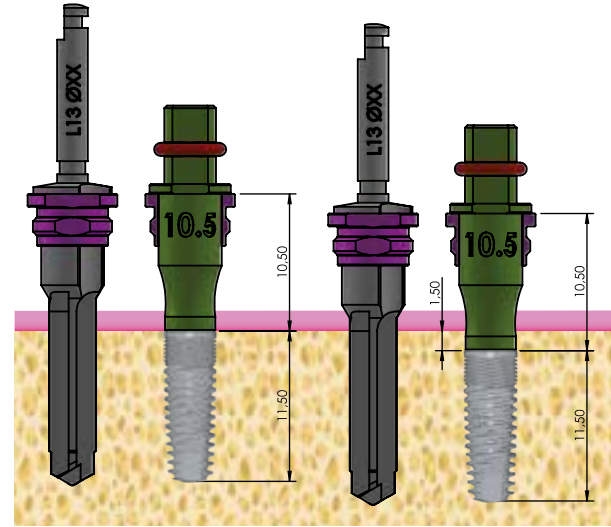


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	10.5	0 - 1.5	13

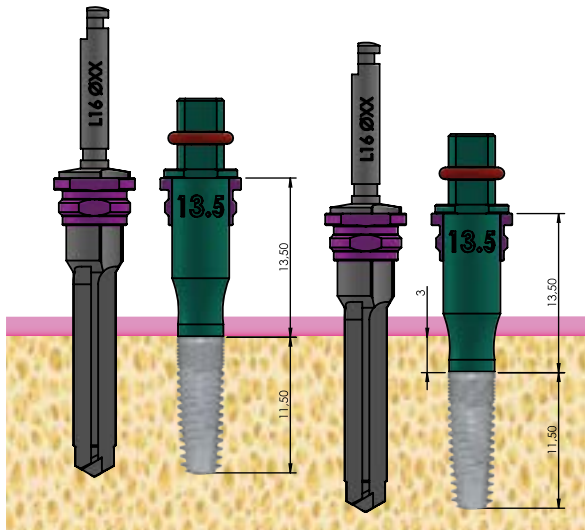


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	13.5	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 13 mm

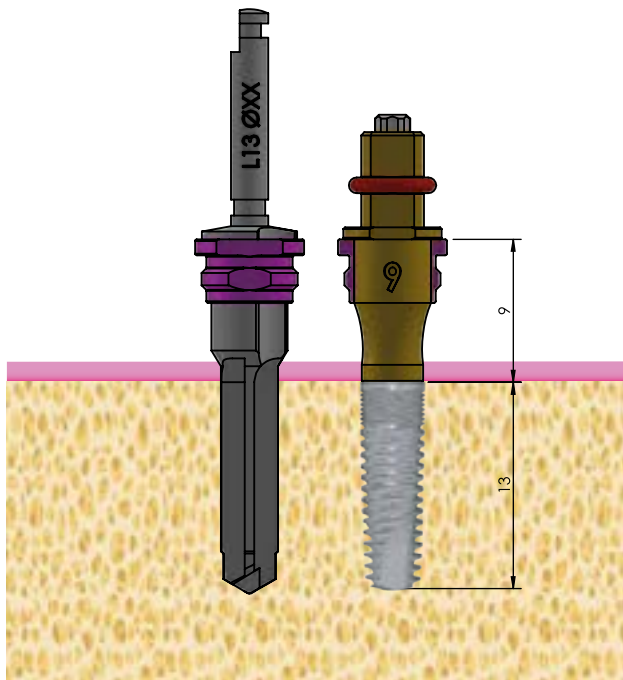


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	9	0	13

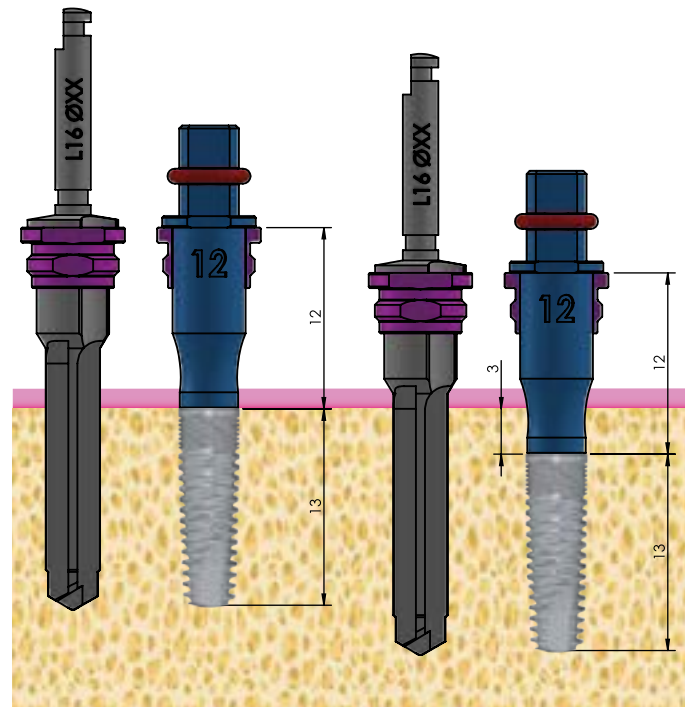


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	12	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 16 mm

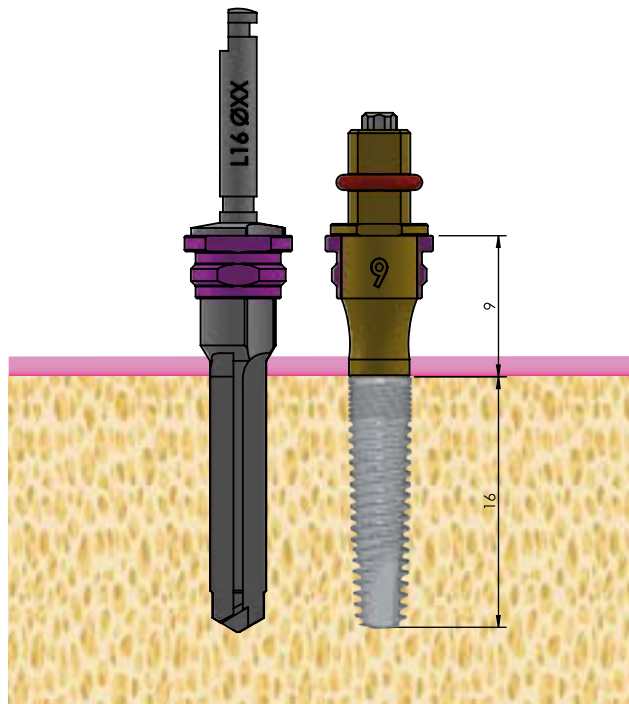

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
16	9	0	16

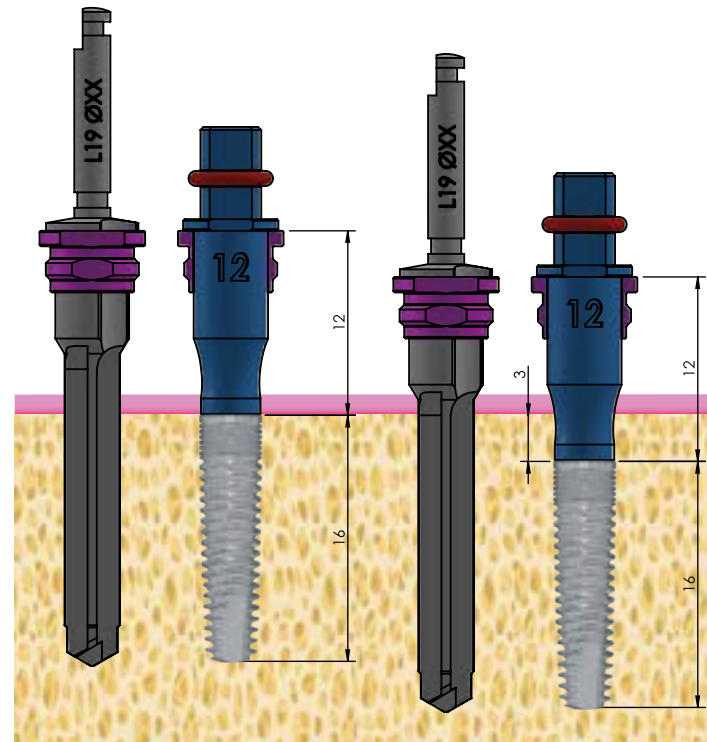

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
16	12	0 - 3	19

Ø XX = Ø Drill


DIO - UF II REGULAR

Implant Length (mm)	Mount (mm)	Depth (mm)	Drill Length (mm)
7	10.5	0 - 1.5	8.5
	12	0 - 3	10
	13.5	0 - 3	11.5
8.5	9	0	8,5
	10.5	0 - 1.5	10
	12	0 - 3	11.5
	13.5	0 - 3	13
10	9	0	10
	10.5	0 - 1.5	11.5
	12	0 - 3	13
11.5	9	0	11.5
	10.5	0 - 1.5	13
	13.5	0 - 3	16
13	9	0	13
	12	0 - 3	16
15	10	0 - 1	16
16	9	0	16
	12	0 - 3	19
18	10	0 - 1	19


DIO		
Model	Ø Implant (mm)	Implant Length (mm)
UF II Regular	3.8	8.5/10/11.5/13/15/16/18
	4/4.5	7/8.5/10/11.5/13/15/16/18

DAS_0030_B1


IMPLANT LENGTH 7 mm

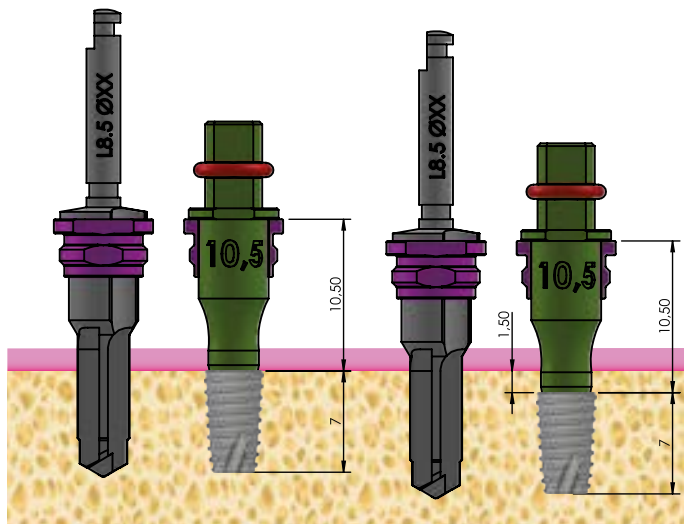

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	10.5	0 - 1.5	8.5

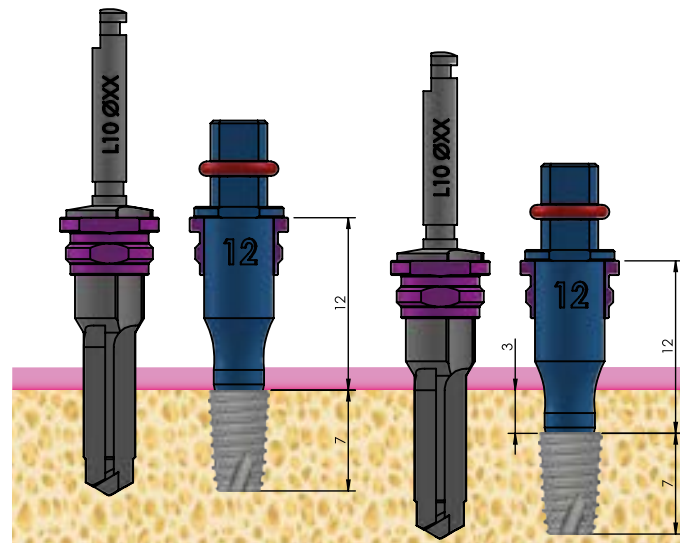


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	12	0 - 3	10

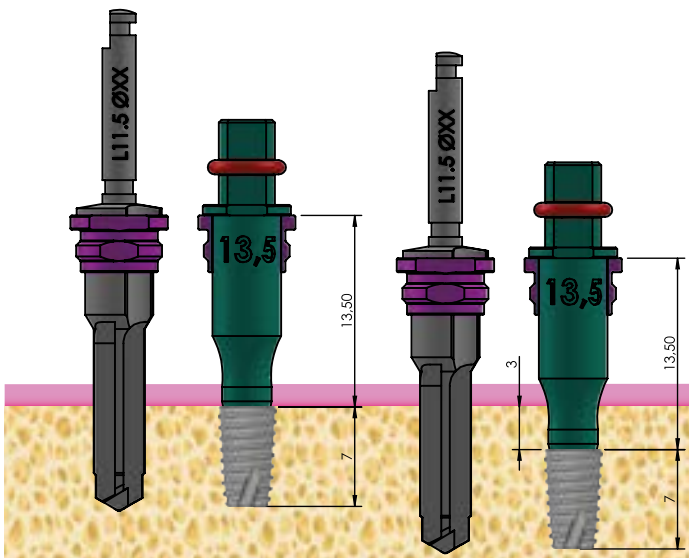


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	13.5	0 - 3	11.5

Ø XX = Ø Drill


IMPLANT LENGTH 8.5 mm

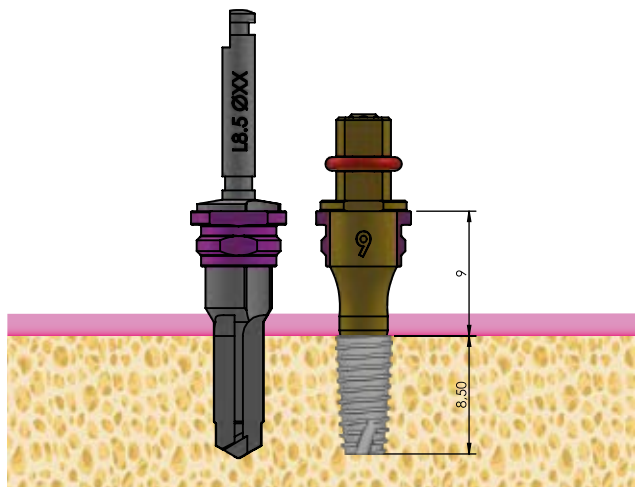


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	9	0	8.5

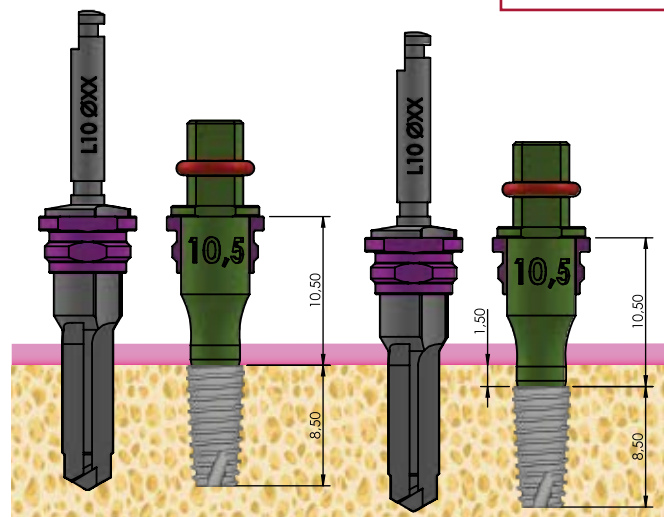 $\emptyset XX = \emptyset \text{ Drill}$ 

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	10.5	0 - 1.5	10

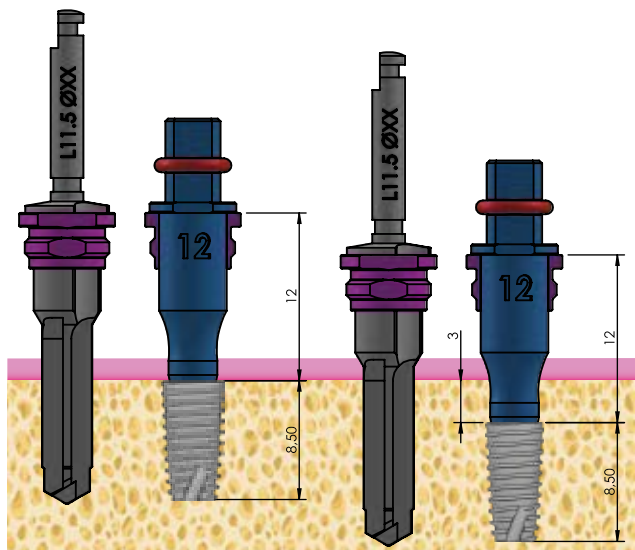


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	12	0 - 3	11.5

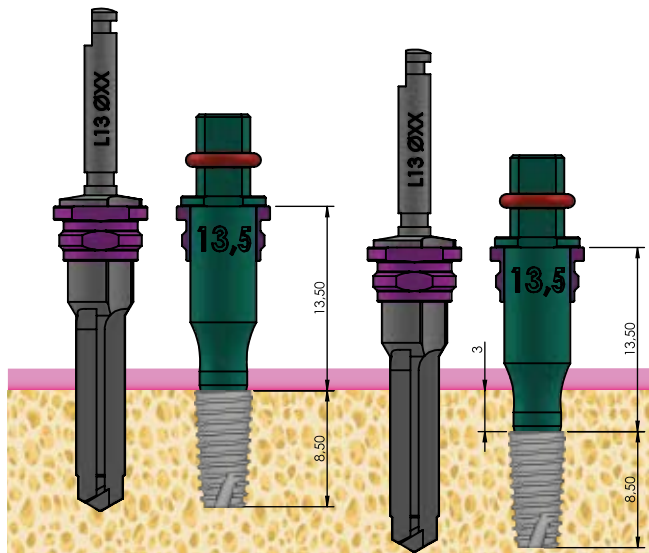


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	13.5	0 - 3	13


IMPLANT LENGTH 10 mm

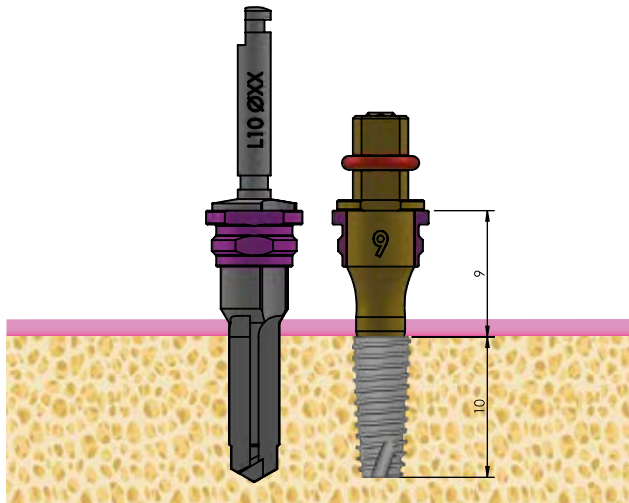

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	9	0	10

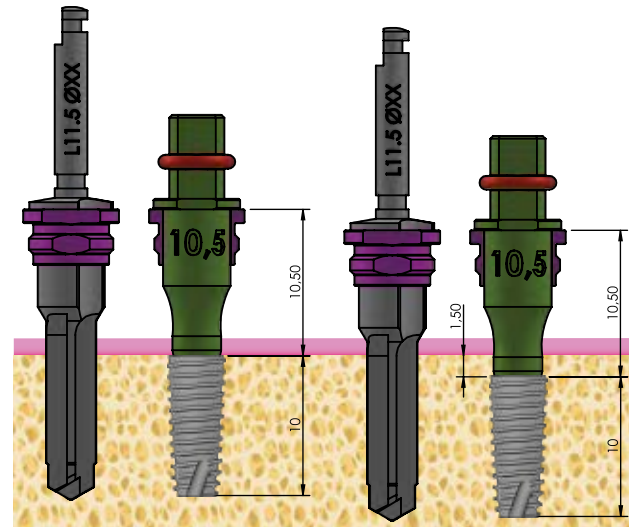


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	10.5	0 - 1.5	11.5

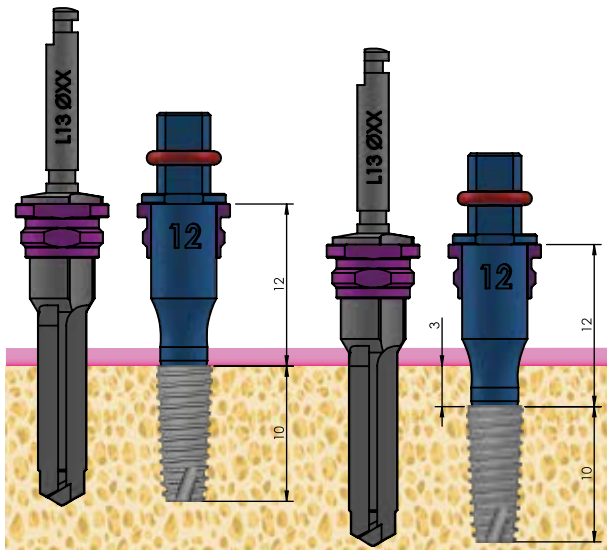


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	12	0 - 3	13

Ø XX = Ø Drill


IMPLANT LENGTH 11.5 mm

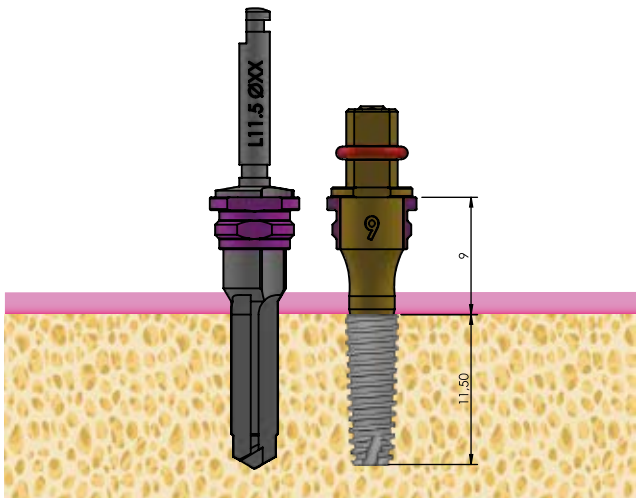

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	9	0	11.5

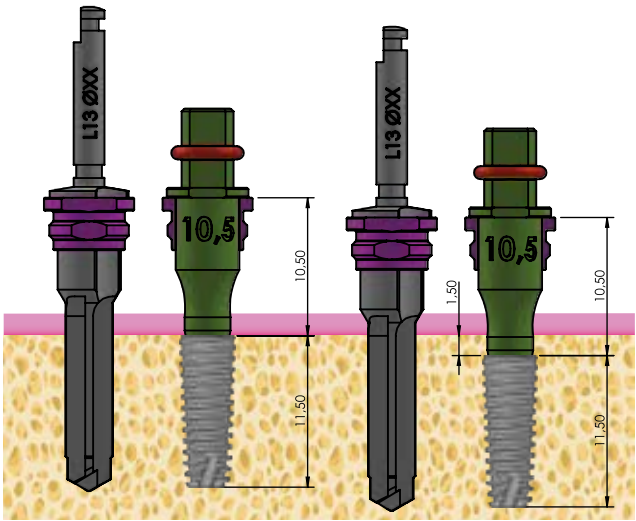


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	10.5	0 - 1.5	13

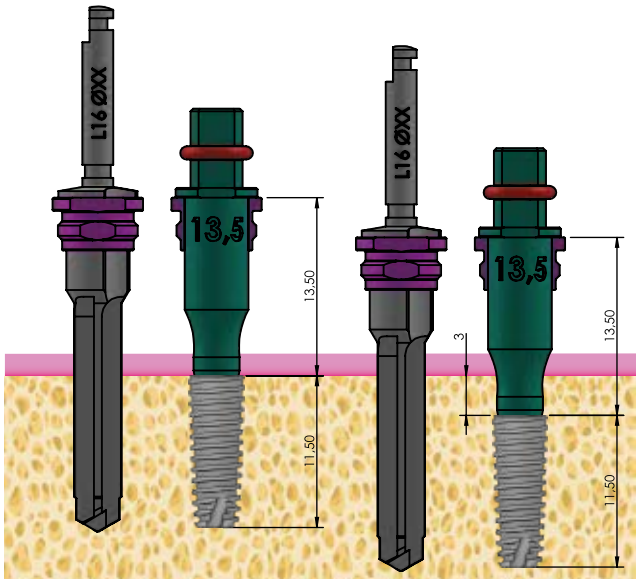


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	13.5	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 13 mm

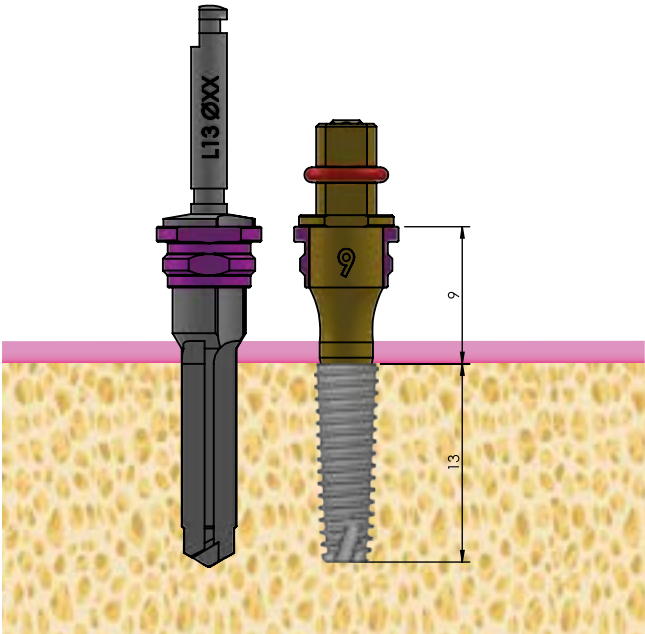


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	9	0	13

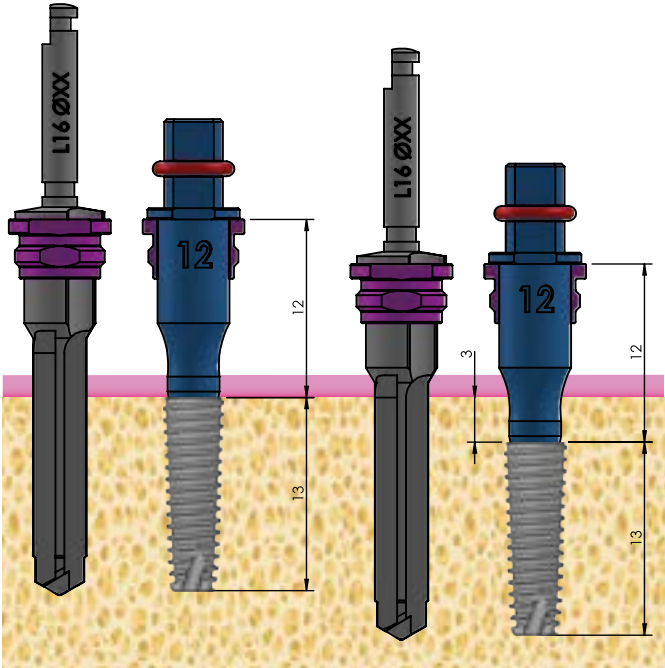

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	12	0 - 3	16


Ø XX = Ø Drill

IMPLANT LENGTH 15 mm

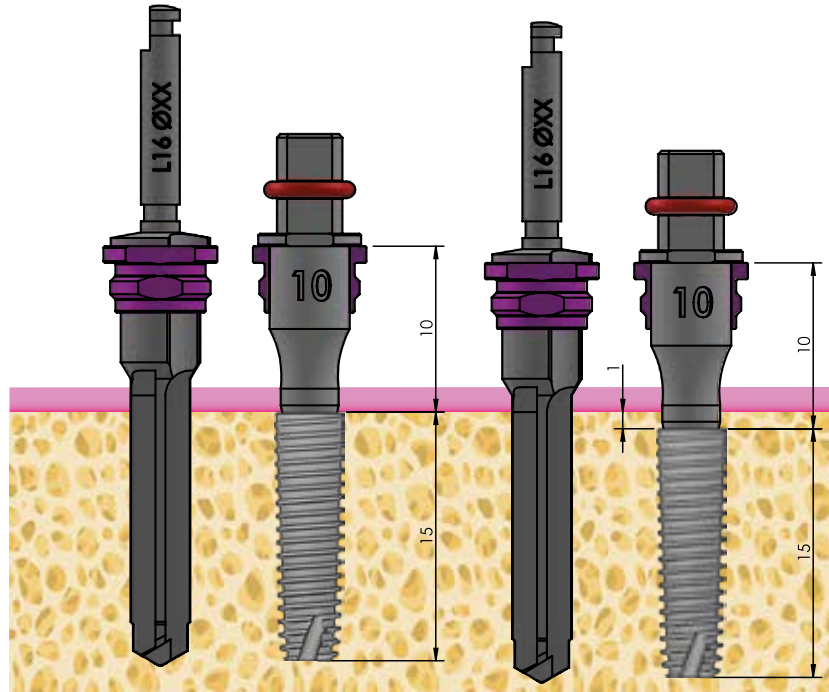


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
15	10	0 - 1	16

IMPLANT LENGTH 16 mm

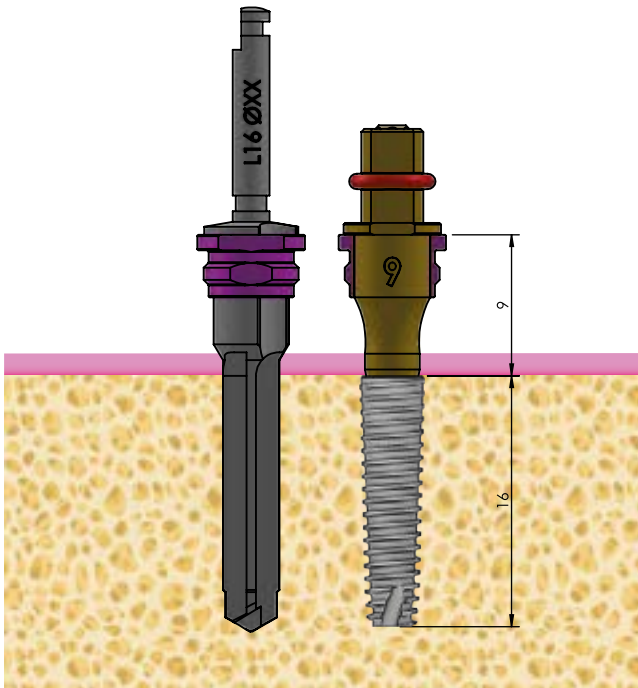


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
16	9	0	16

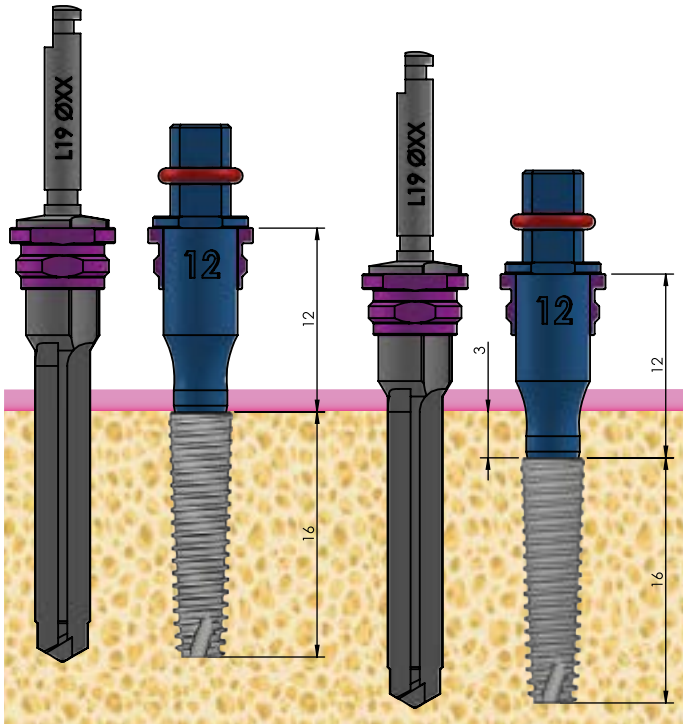

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
16	12	0 - 3	19

Ø XX = Ø Drill

IMPLANT LENGTH 18 mm

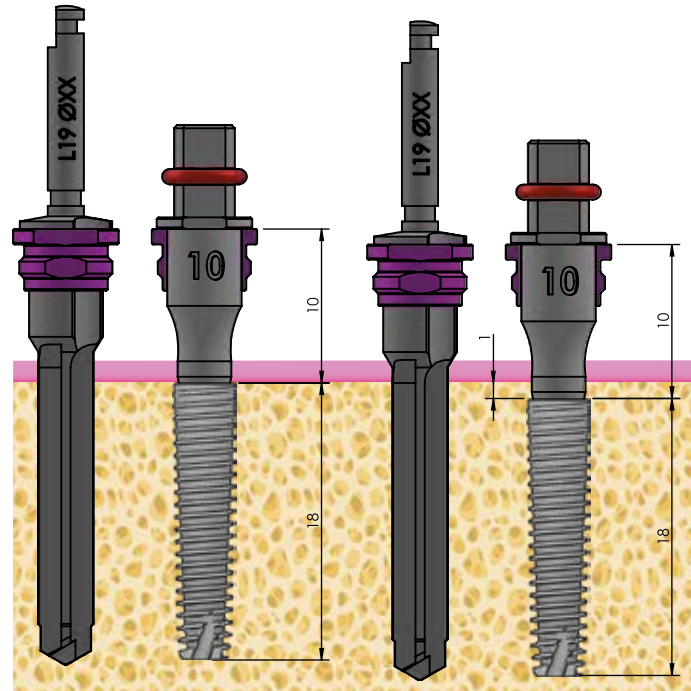

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
18	10	0 - 1	19

Ø XX = Ø Drill

NOTES


HIOSSEN - ET SYSTEM

Implant Length (mm)	Mount (mm)	Depth (mm)	Drill Length (mm)
7	10.5	0 - 1.5	8.5
	12	0 - 3	10
	13.5	0 - 3	11.5
8.5	9	0	8.5
	10.5	0 - 1.5	10
	12	0 - 3	11.5
	13.5	0 - 3	13
10	9	0	10
	10.5	0 - 1.5	11.5
	12	0 - 3	13
11.5	9	0	11.5
	10.5	0 - 1.5	13
	13.5	0 - 3	16
13	9	0	13
	12	0 - 3	16
15	10	0 - 1	16


HIOSSEN ET SYSTEM	
Ø Implant (mm)	Implant Length (mm)
4/4.5	7/8.5/10/11.5/13/15

DAS_0030_B1


IMPLANT LENGTH 7 mm

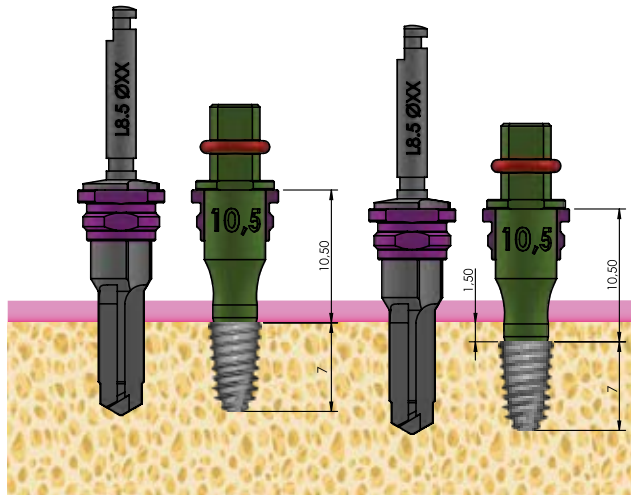

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	10.5	0 - 1.5	8.5

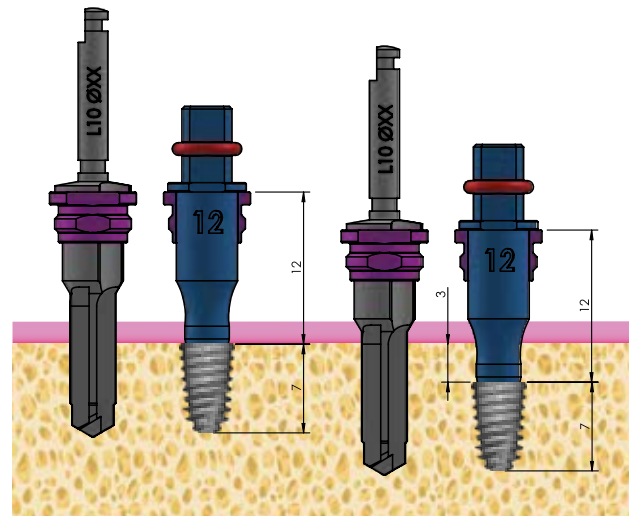


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	12	0 - 3	10

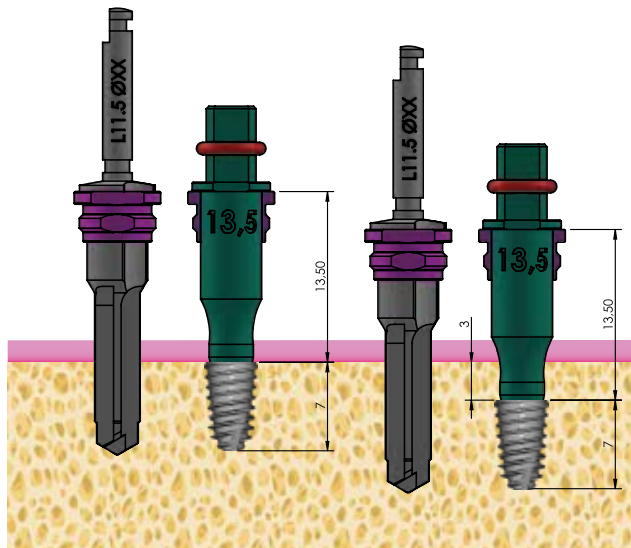


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	13.5	0 - 3	11.5

Ø XX = Ø Drill


IMPLANT LENGTH 8.5 mm

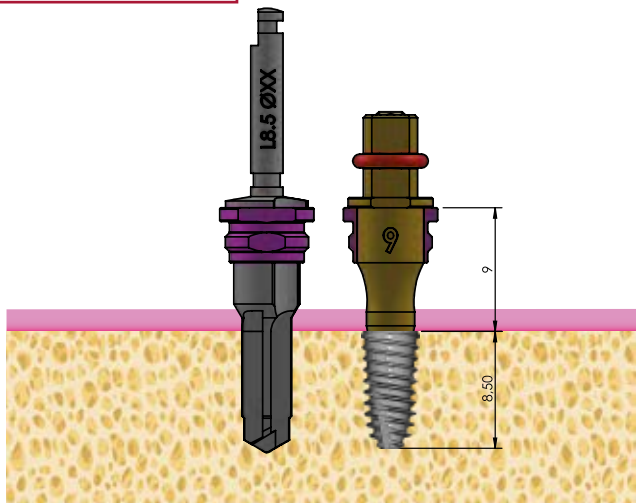


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	9	0	8.5

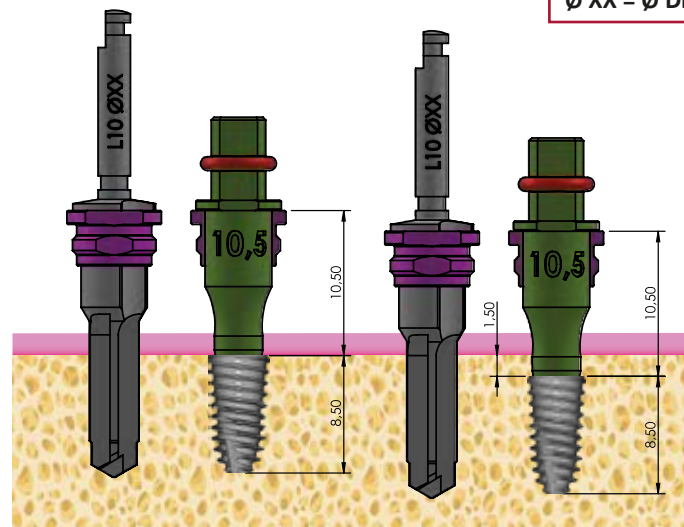 $\varnothing XX = \varnothing \text{ Drill}$


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	10.5	0 - 1.5	10

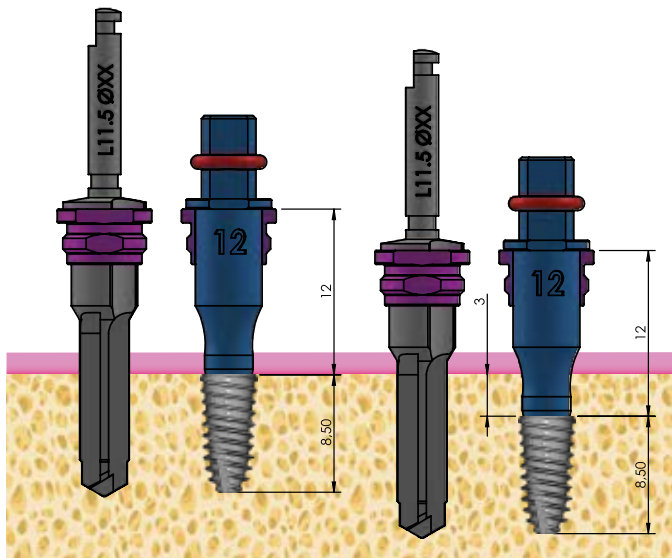


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	12	0 - 3	11.5

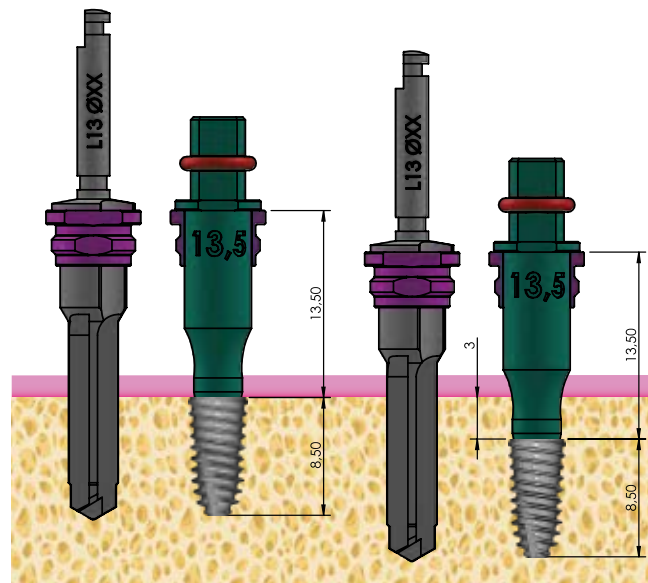


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	13.5	0 - 3	13


IMPLANT LENGTH 10 mm

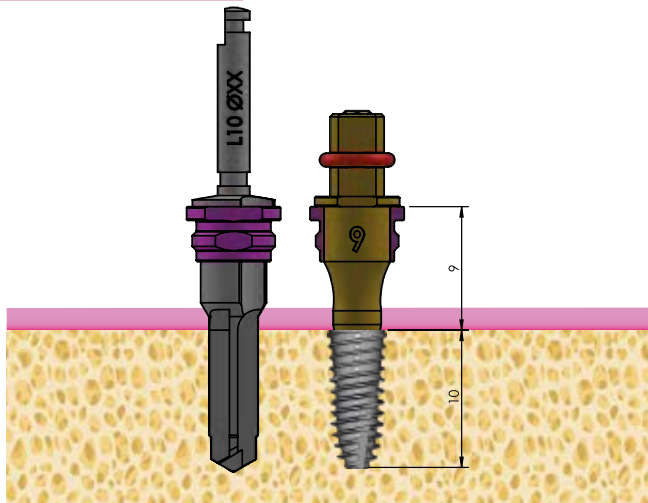

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	9	0	10

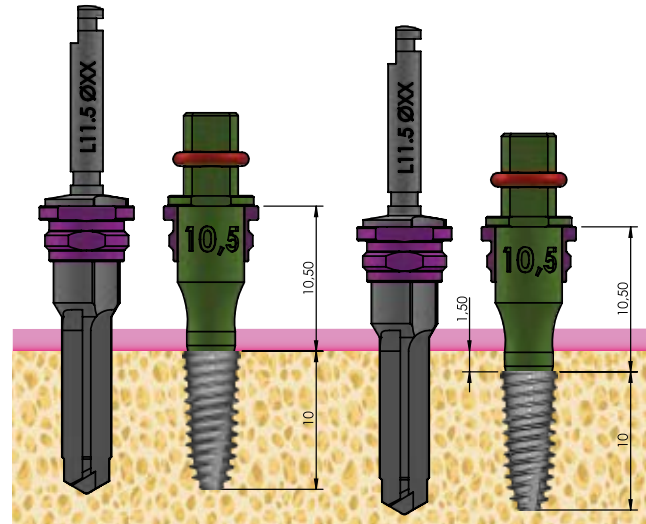


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	10.5	0 - 1.5	11.5

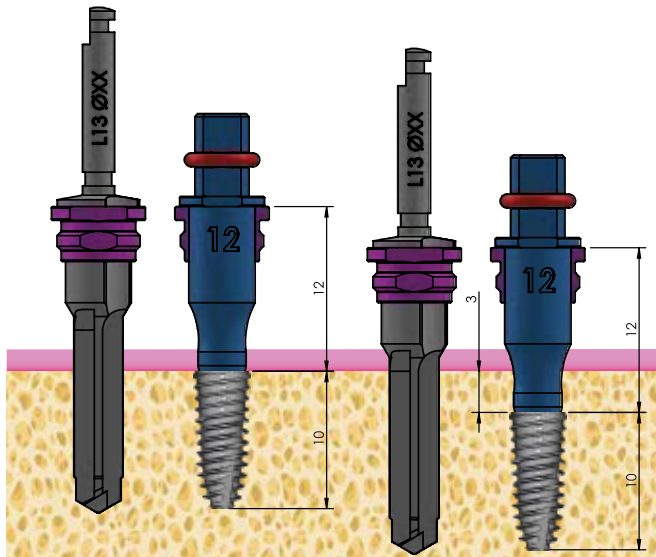


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	12	0 - 3	13

Ø XX = Ø Drill


IMPLANT LENGTH 11.5 mm

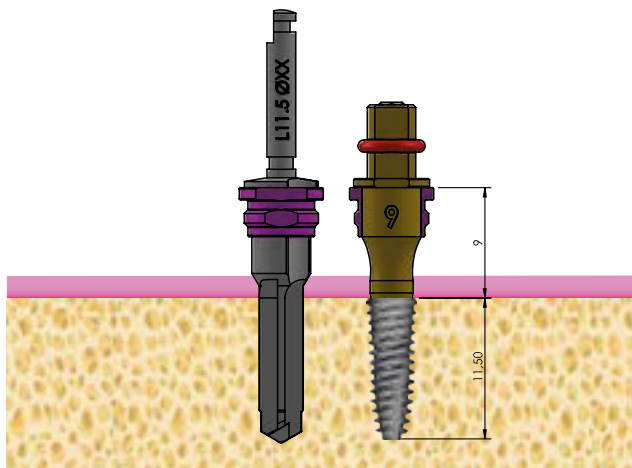

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	9	0	11.5

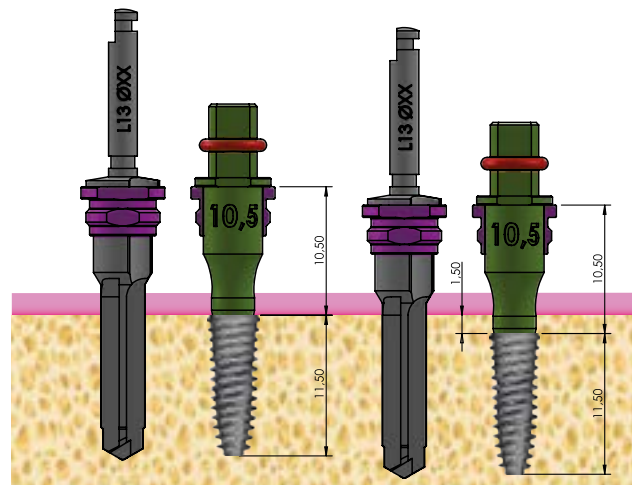


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	10.5	0 - 1.5	13

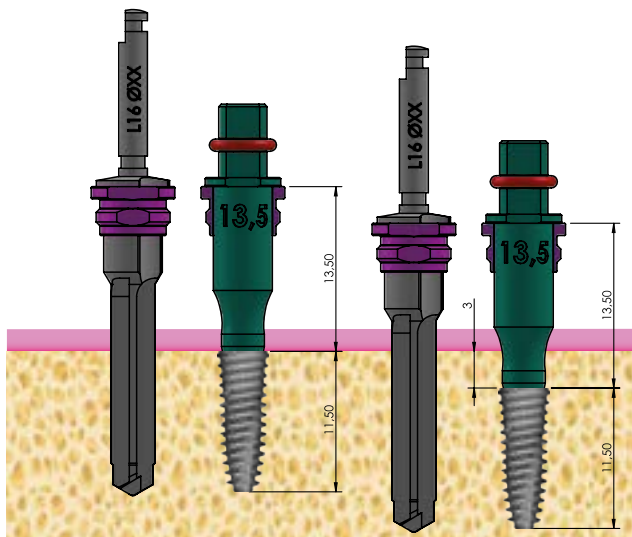


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	13.5	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 13 mm

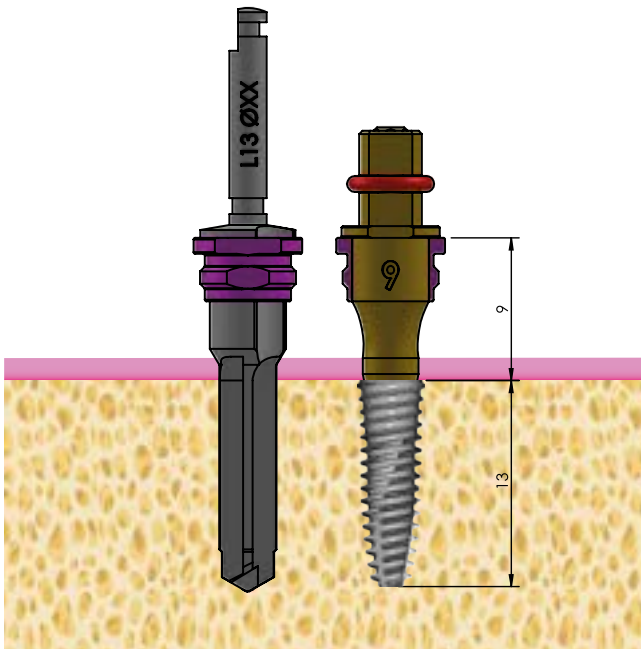


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	9	0	13

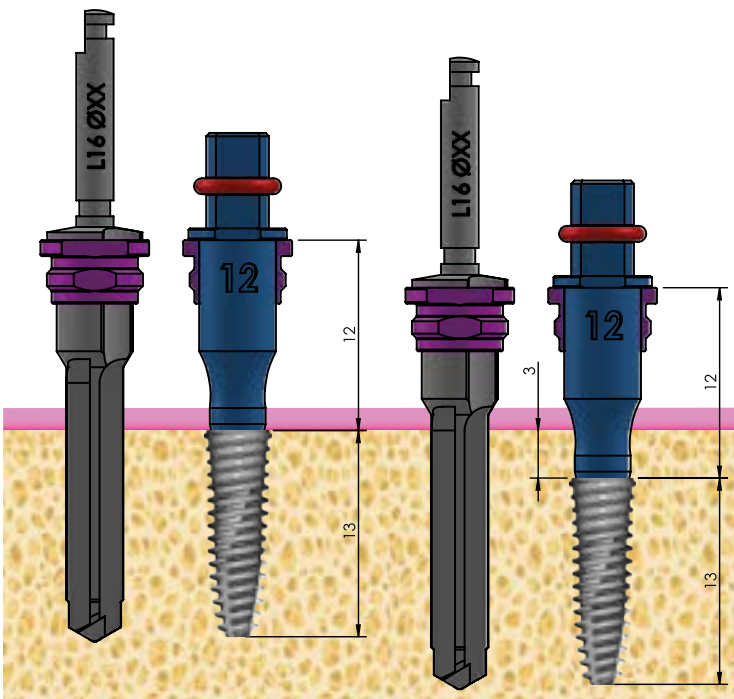

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	12	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 15 mm

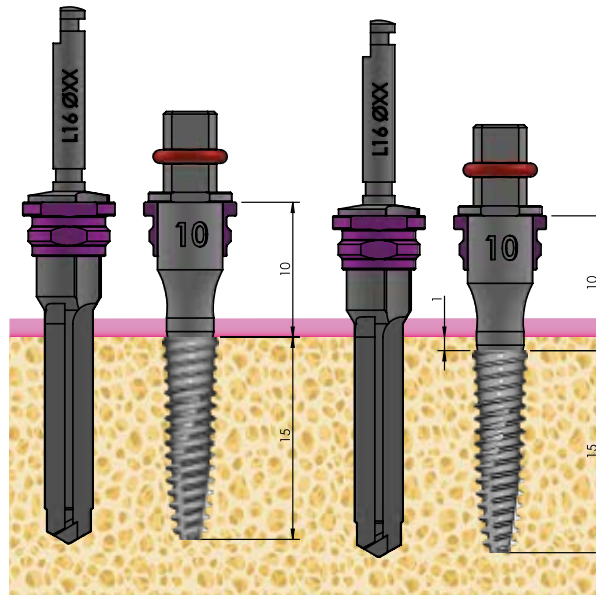

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
15	10	0 - 1	16

Ø XX = Ø Drill


NOTES

MIS INTERNAL HEXAGON - SEVEN/M4

Implant Length (mm)	Mount (mm)	Depth (mm)	Drill Length (mm)
6	9	0	6
8	9.5	0 - 0.5	8.5
10	9	0	10
	10.5	0 - 1.5	11.5
	12	0 - 3	13
11.5	9	0	11.5
	10.5	0 - 1.5	13
	13.5	0 - 3	16
13	9	0	13
	12	0 - 3	16
16	9	0	16
	12	0 - 3	19


MIS Internal Hexagon		
Model	Ø Implant (mm)	Implant Length (mm)
Seven SP	3.75	8/10/11.5/13/16
Seven SP	4.2	6/8/10/11.5/13/16
M4 SP	3.75	8/10/11.5/13/16
M4 SP	4.2	6/8/10/11.5/13/16

IMPLANT LENGTH 6 mm

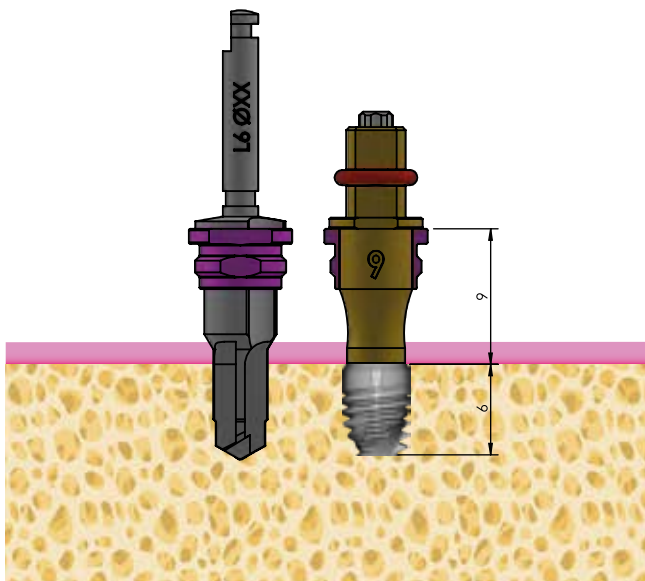

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
6	9	0	6


IMPLANT LENGTH 8 mm

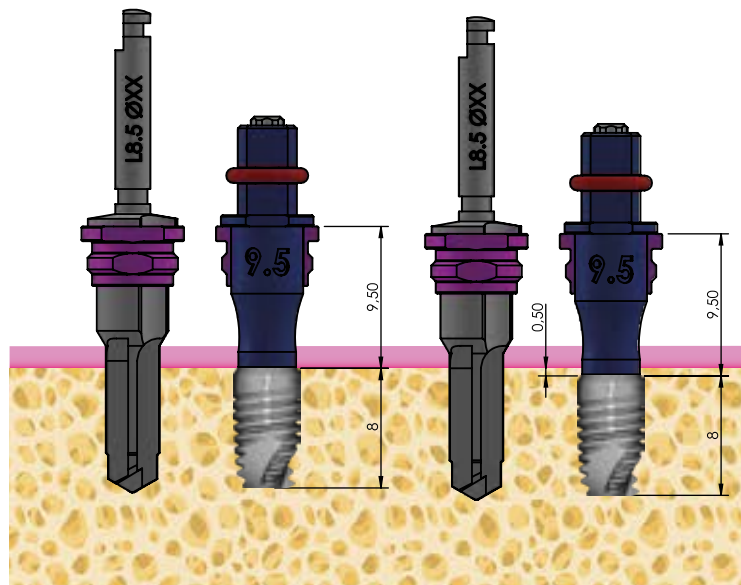


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8	9.5	0 - 0.5	8.5

Ø XX = Ø Drill


IMPLANT LENGTH 10 mm

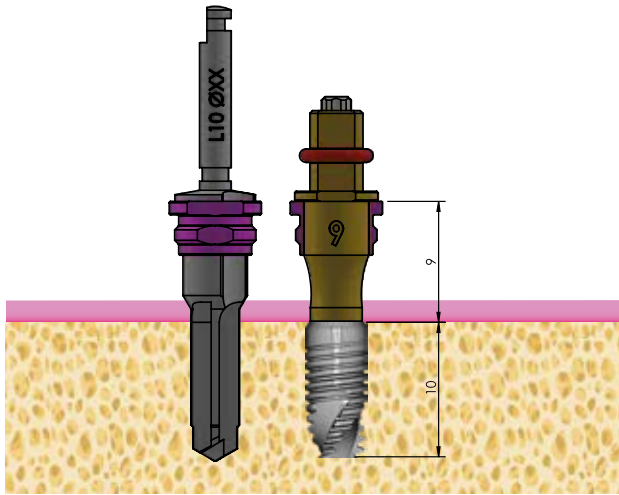

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	9	0	10

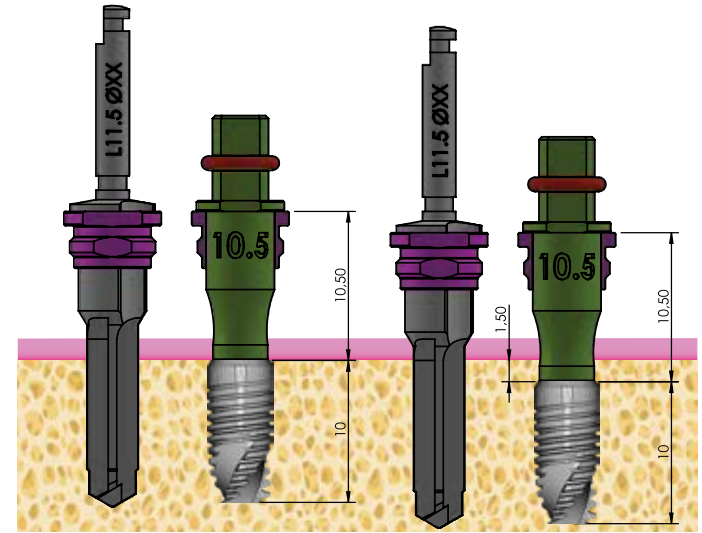


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	10.5	0 - 1.5	11.5

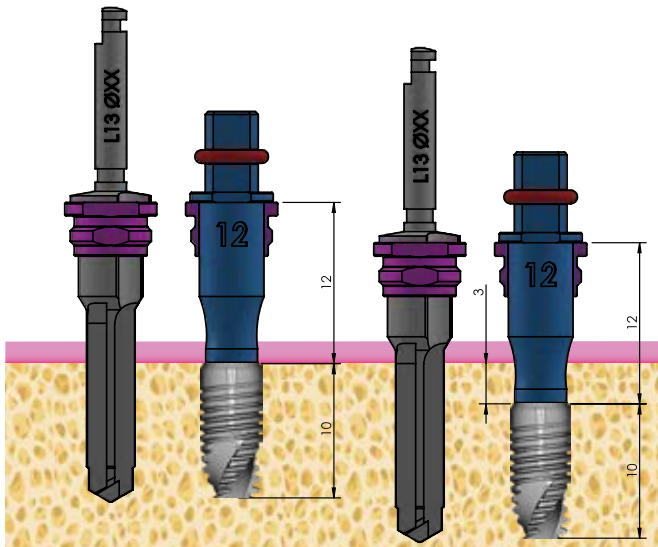


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	12	0 - 3	13

 $\emptyset XX = \emptyset \text{ Drill}$


IMPLANT LENGTH 11.5 mm

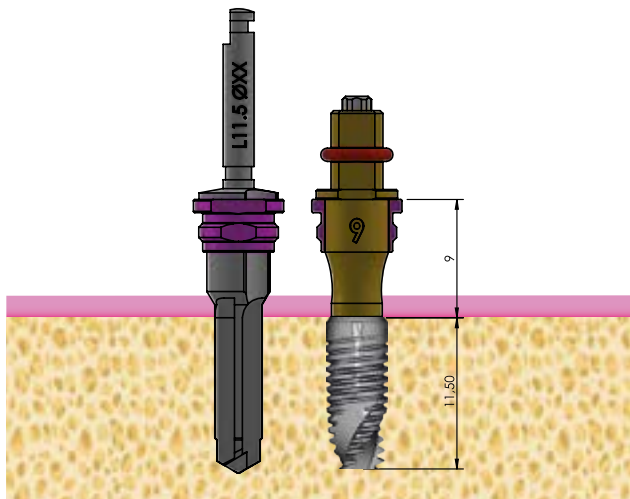

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	9	0	11.5

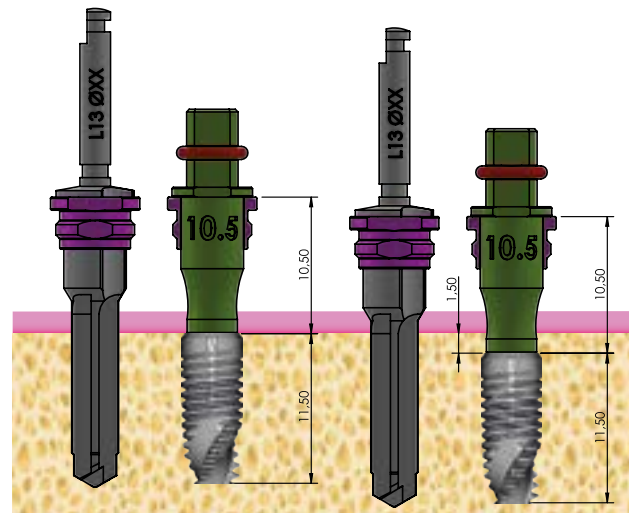


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	10.5	0 - 1.50	13

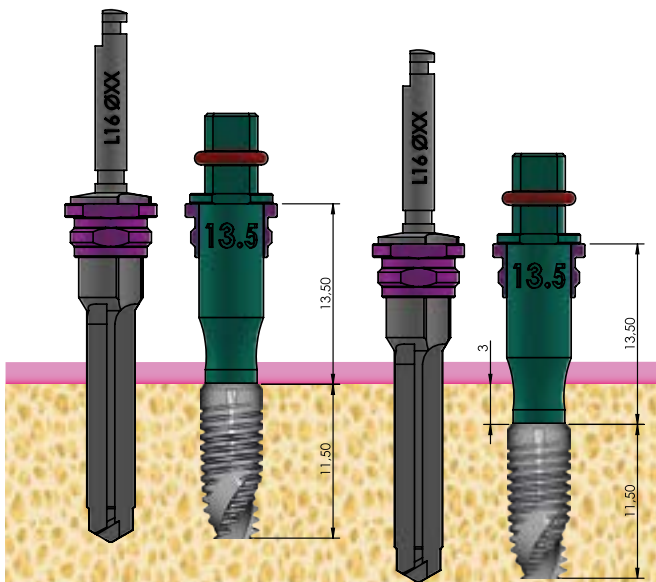


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	13.5	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 13 mm

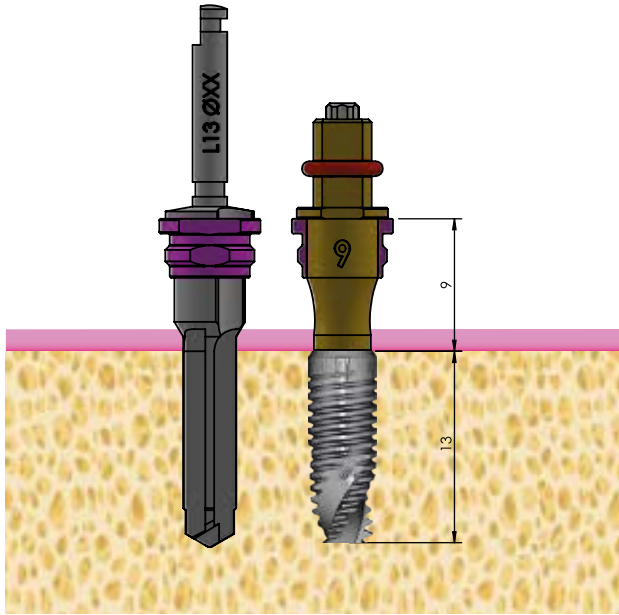


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	9	0	13

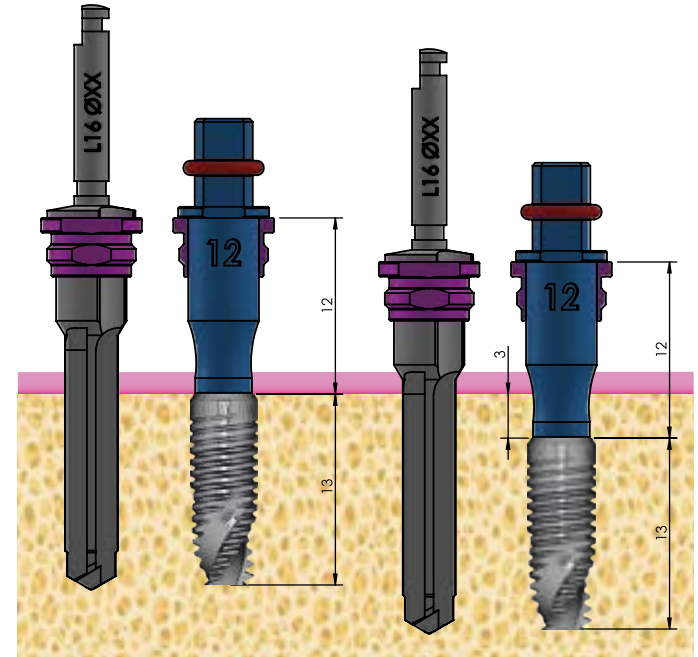


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	12	0 - 3	16

 $\varnothing XX = \varnothing \text{ Drill}$

IMPLANT LENGTH 16 mm

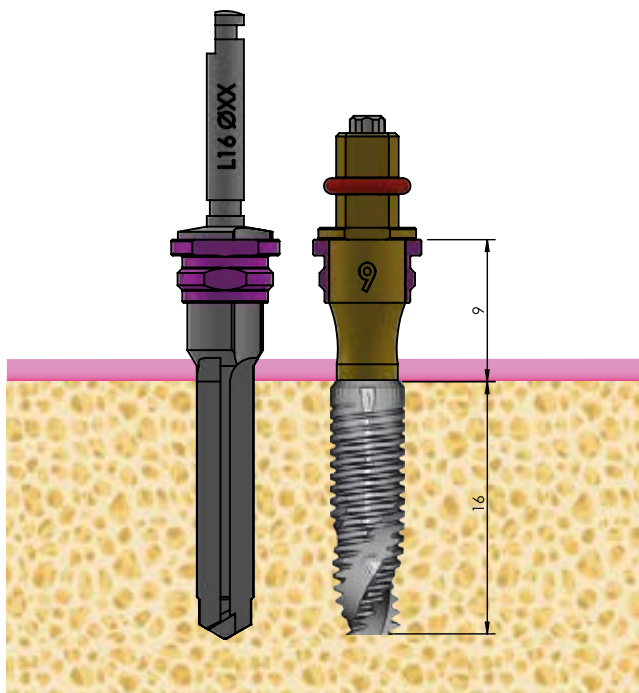

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
16	9	0	16

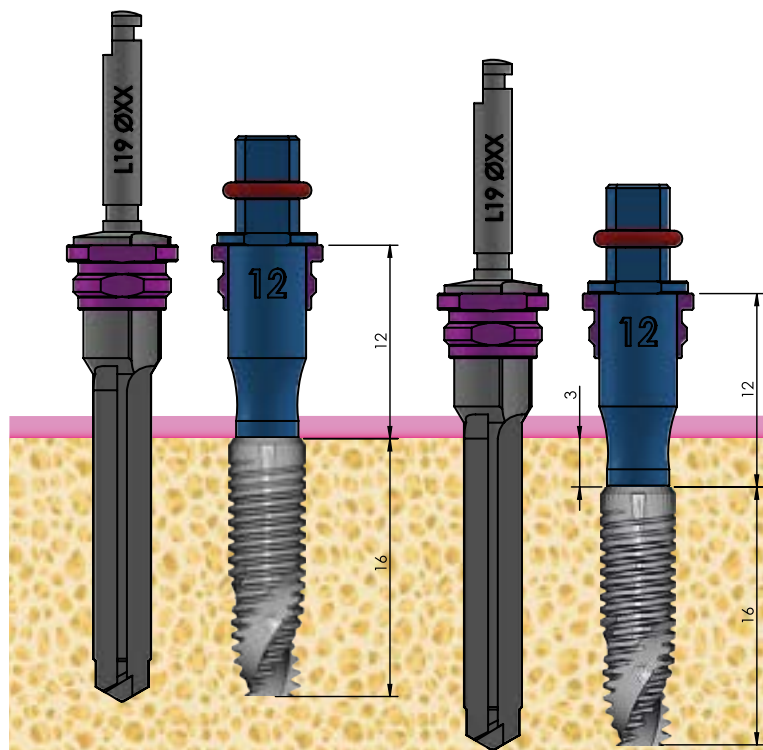

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
16	12	0 - 3	19

Ø XX = Ø Drill


NEOBIOTECH - IS SYSTEM

Implant Length (mm)	Mount (mm)	Depth (mm)	Drill Length (mm)
8.5	9	0	8.5
	10.5	0 - 1.5	10
	12	0 - 3	11.5
	13.5	0 - 3	13
10	9	0	10
	10.5	0 - 1.5	11.5
	12	0 - 3	13
11.5	9	0	11.5
	10.5	0 - 1.5	13
	13.5	0 - 3	16
13	9	0	13
	12	0 - 3	16

NEOBIOTECH IS SYSTEM	
Ø Implant (mm)	Implant Length (mm)
3.5/4/4.5	8.5/10/11.5/13


DAS_0030_B1

IMPLANT LENGTH 8.5 mm

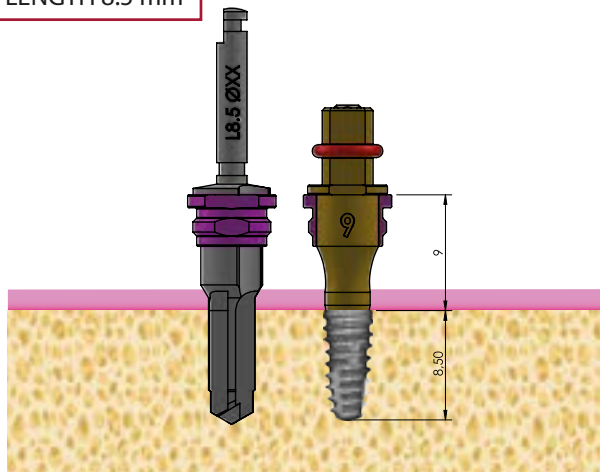


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	9	0	8.5


Ø XX = Ø Drill

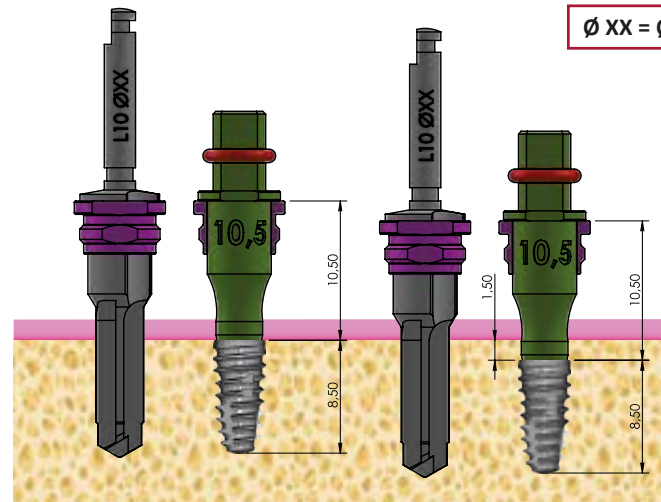

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	10.5	0 - 1.5	10

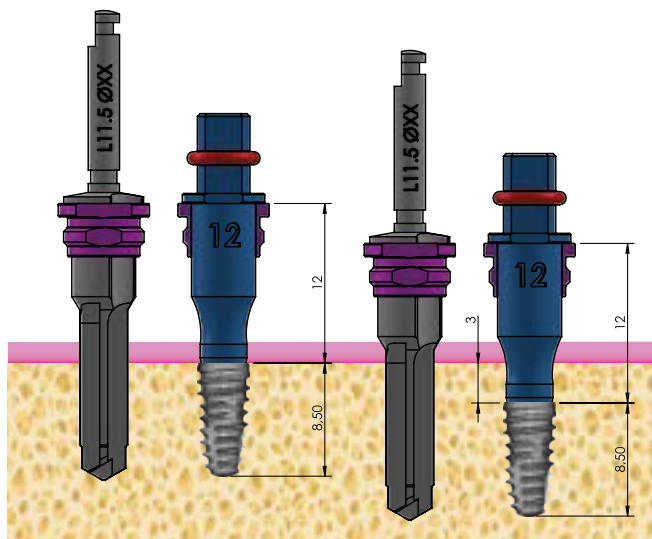


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	12	0 - 3	11.5

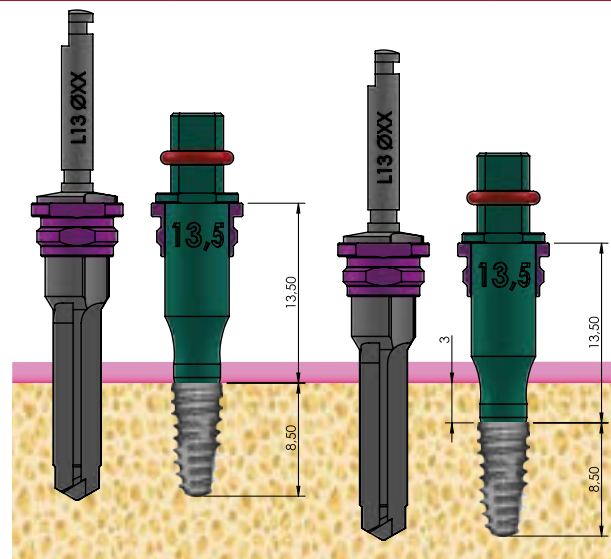


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	13.5	0 - 3	13


IMPLANT LENGTH 10 mm

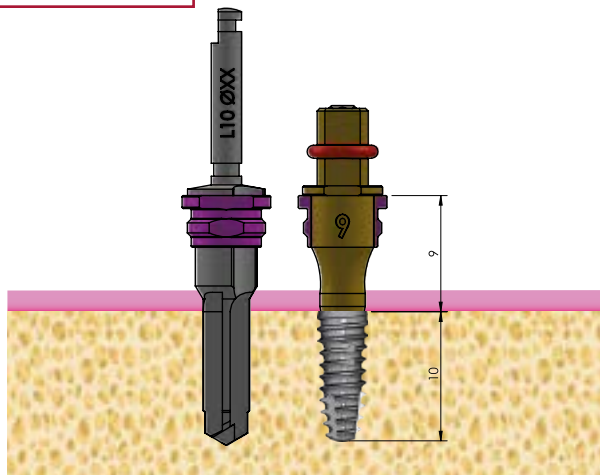

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	9	0	10

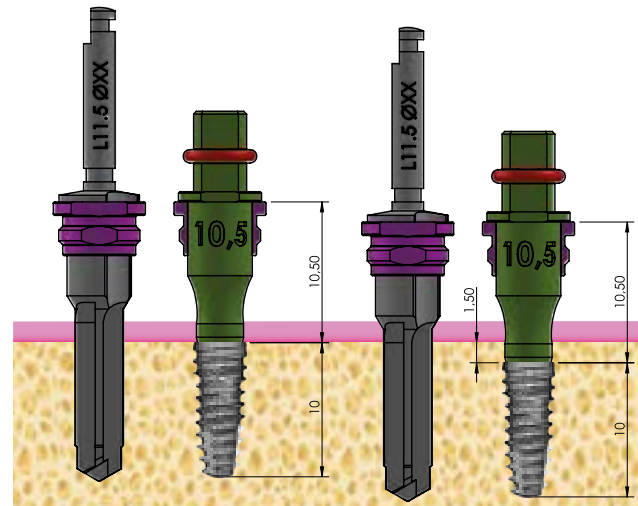


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	10.5	0 - 1.5	11.5

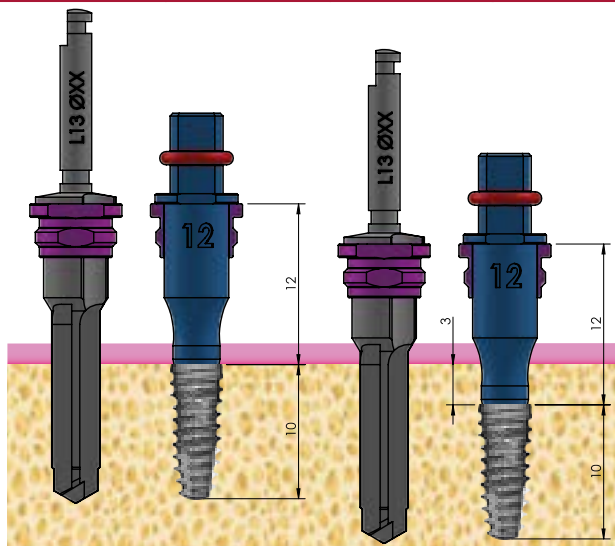


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	12	0 - 3	13

Ø XX = Ø Drill


IMPLANT LENGTH 11.5 mm

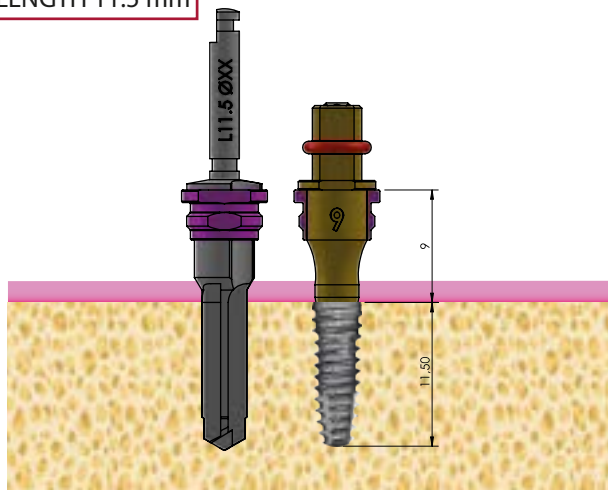

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	9	0	11.5

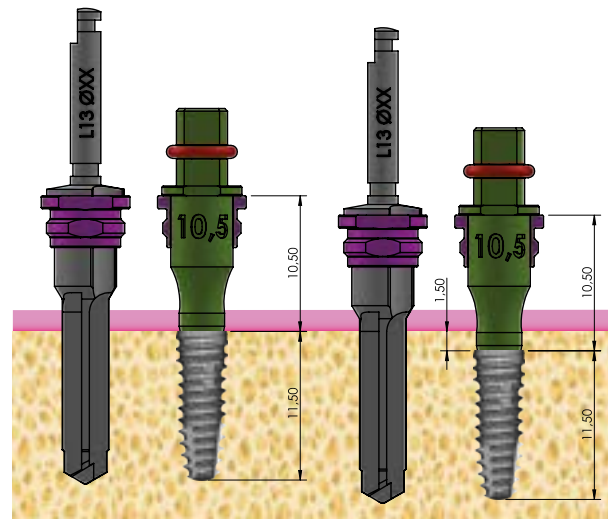


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	10.5	0 - 1.5	13

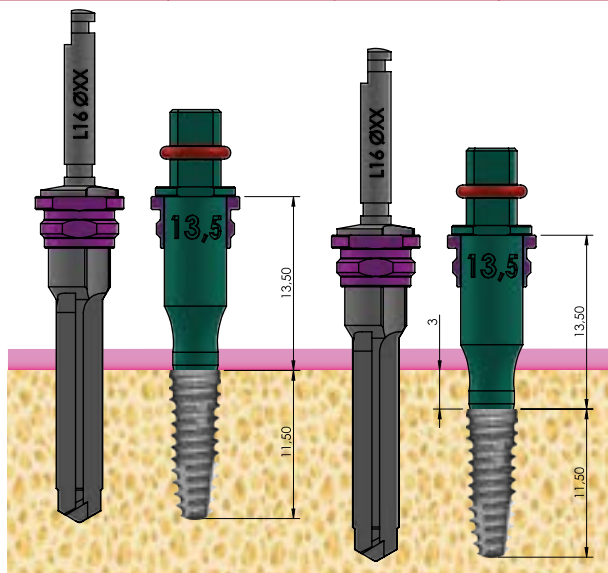


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	13.5	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 13 mm

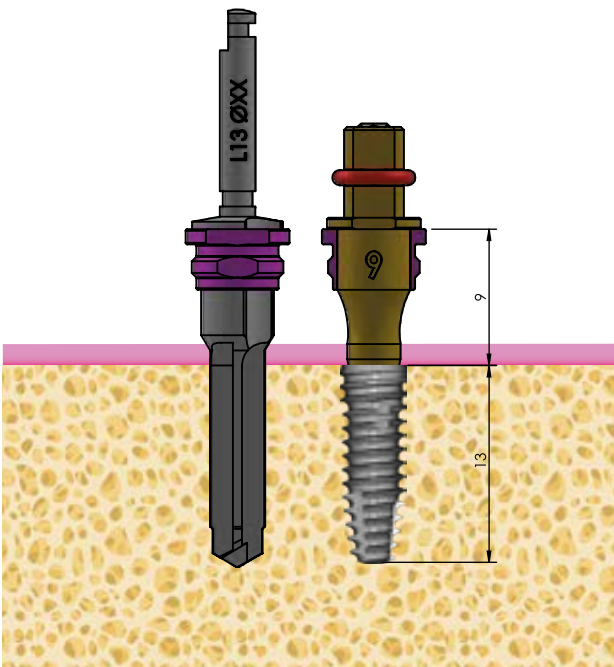

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	9	0	13

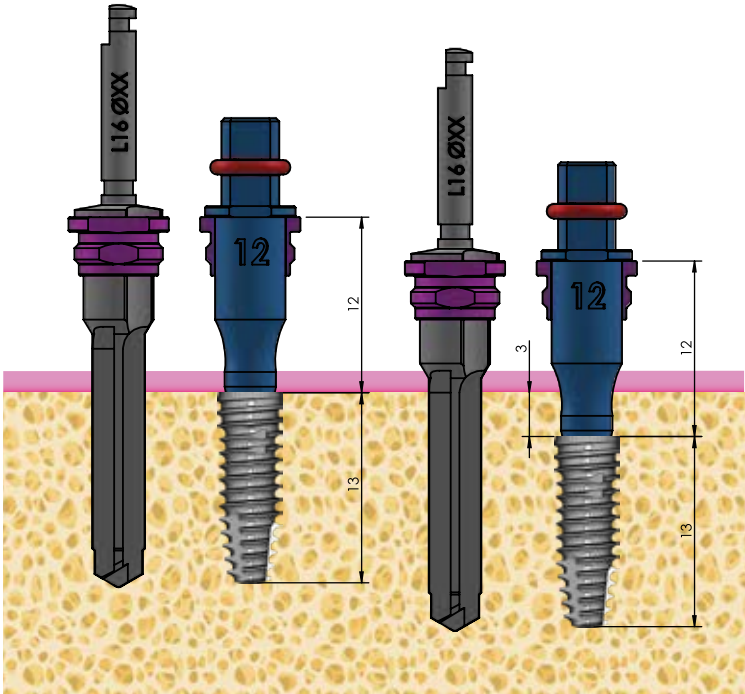

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	12	0 - 3	16

Ø XX = Ø Drill

NOTES

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.


NOBEL BIO CARE - ACTIVE/PARALLEL CONICAL CONNECTION/REPLACE CONICAL CONNECTION

Implant Length (mm)	Mount (mm)	Depth (mm)	Drill Length (mm)
7	10.5	0 - 1.5	8.5
	12	0 - 3	10
	13.5	0 - 3	11.5
8.5	9	0	8.5
	10.5	0 - 1.5	10
	12	0 - 3	11.5
	13.5	0 - 3	13
10	9	0	10
	10.5	0 - 1.5	11.5
	12	0 - 3	13
11.5	9	0	11.5
	10.5	0 - 1.5	13
	13.5	0 - 3	16
13	9	0	13
	12	0 - 3	16
15	10	0 - 1	16
16	9	0	16
	12	0 - 3	19
18	10	0 - 1	19


NOBEL		
Model	Ø Implant (mm)	Implant Length (mm)
Active	3.5	8.5/10/11.5/13/15/18
Parallel Conical Connection	3.75	7/8.5/10/11.5/13/15/18
Replace Conical Connection	3.5	10/11.5/13/16

DAS_0021_B1


IMPLANT LENGTH 7 mm

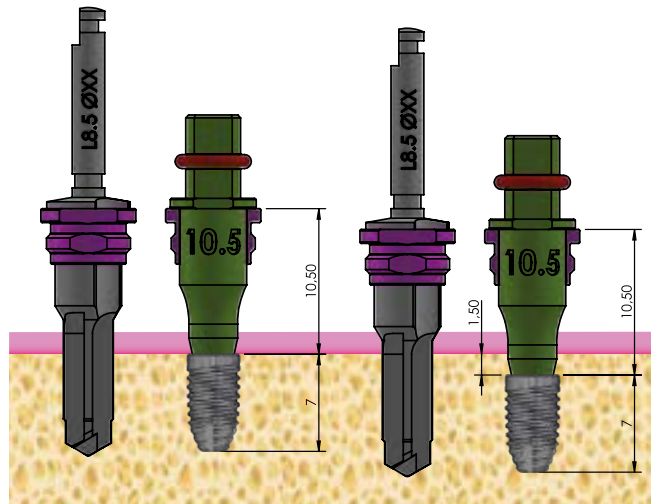

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	10.5	0 - 1.5	8.5

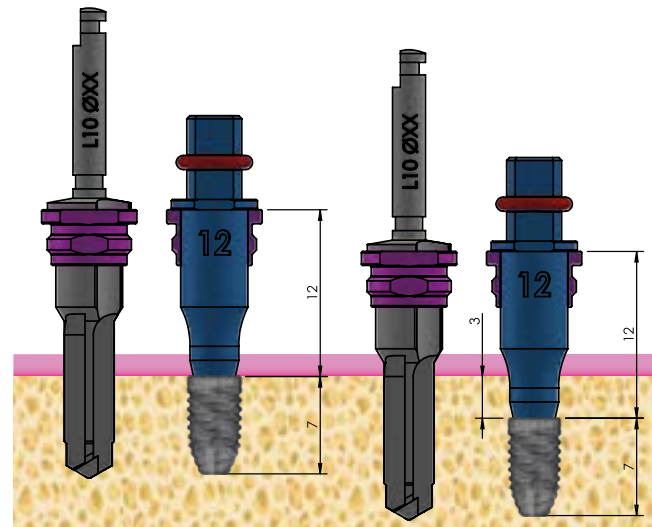


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	12	0 - 3	10

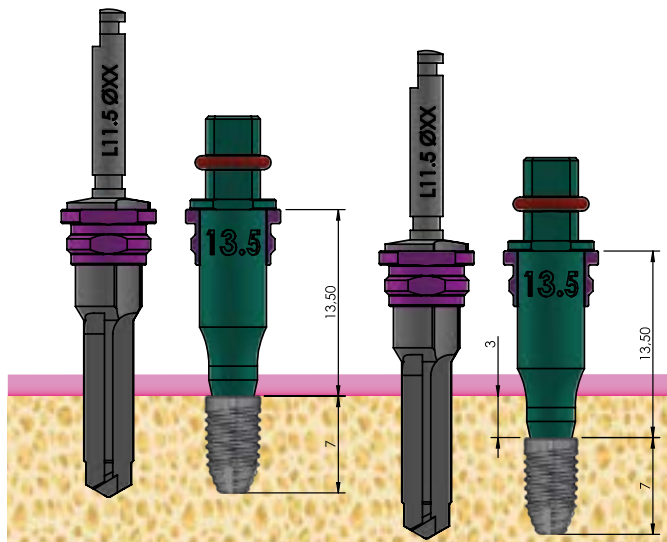

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	13.5	0 - 3	11.5


Ø XX = Ø Drill

IMPLANT LENGTH 8.5 mm

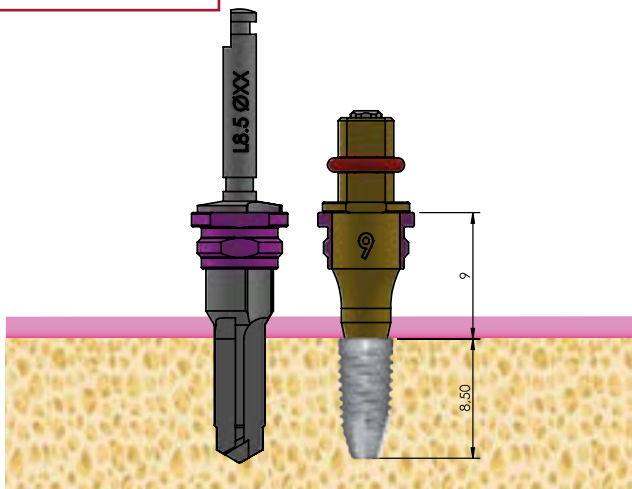


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	9	0	8.5


Ø XX = Ø Drill

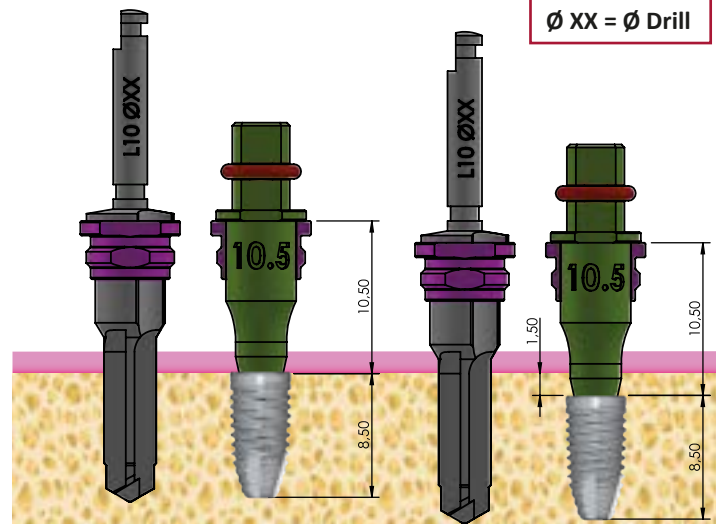

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	10.5	0 - 1.5	10

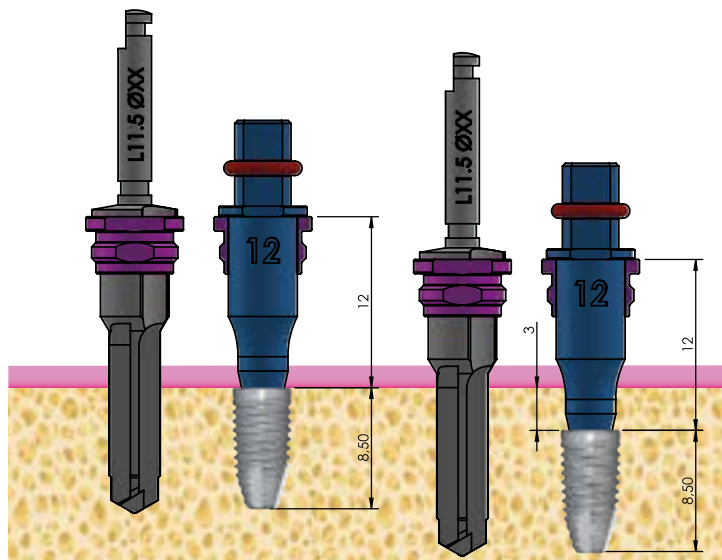


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	12	0 - 3	11.5

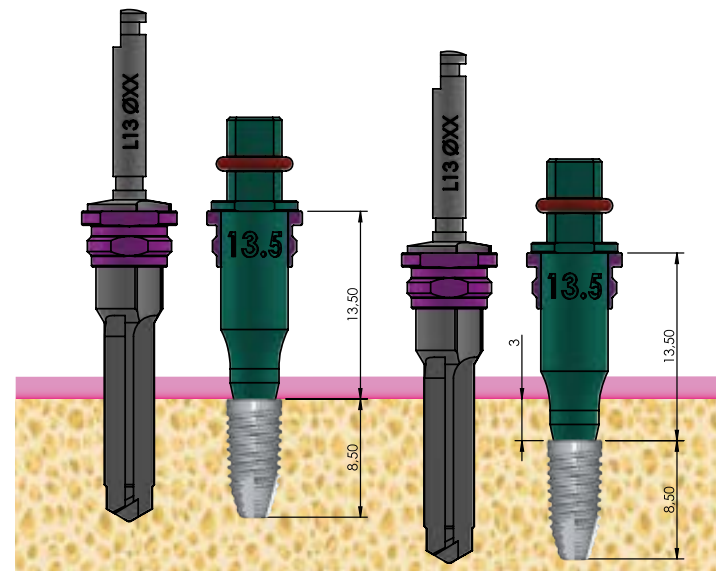


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8.5	13.5	0 - 3	13


IMPLANT LENGTH 10 mm

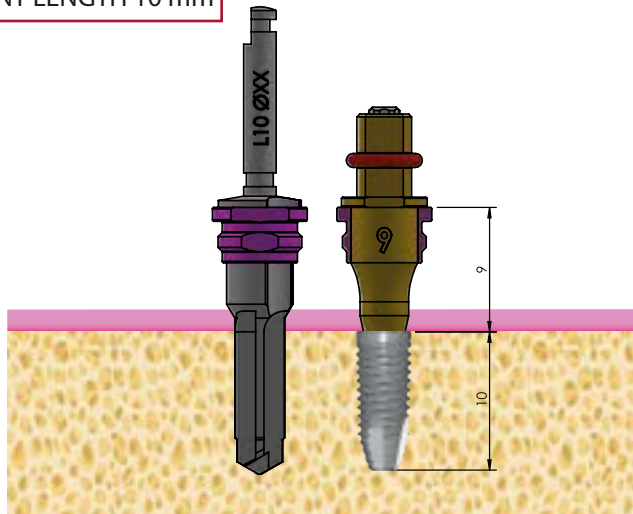

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	9	0	10

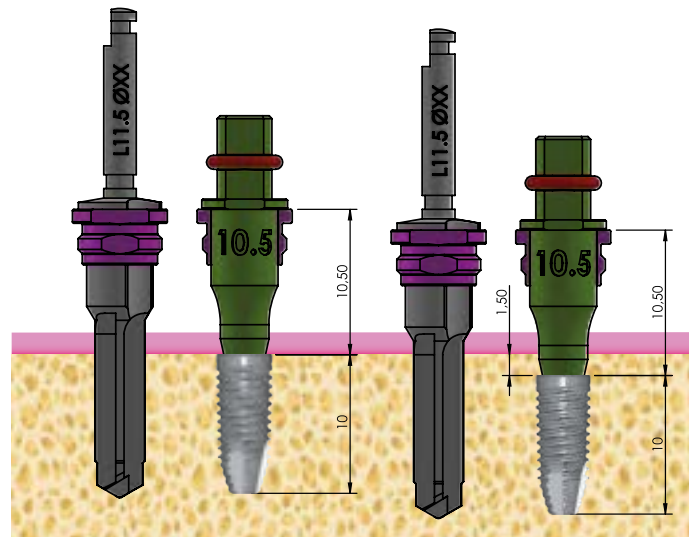


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	10.5	0 - 1.50	11.5

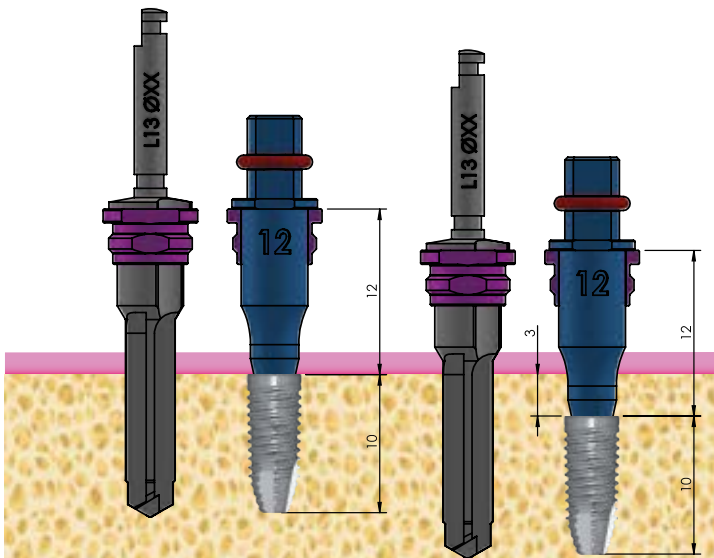


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	12	0 - 3	13

Ø XX = Ø Drill


IMPLANT LENGTH 11.5 mm

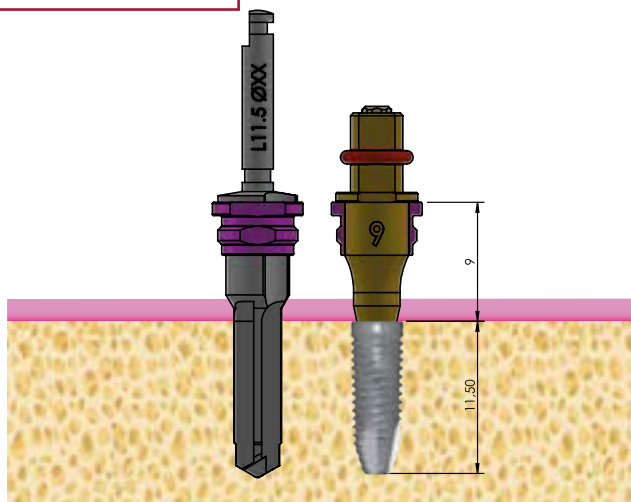

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	9	0	11.5

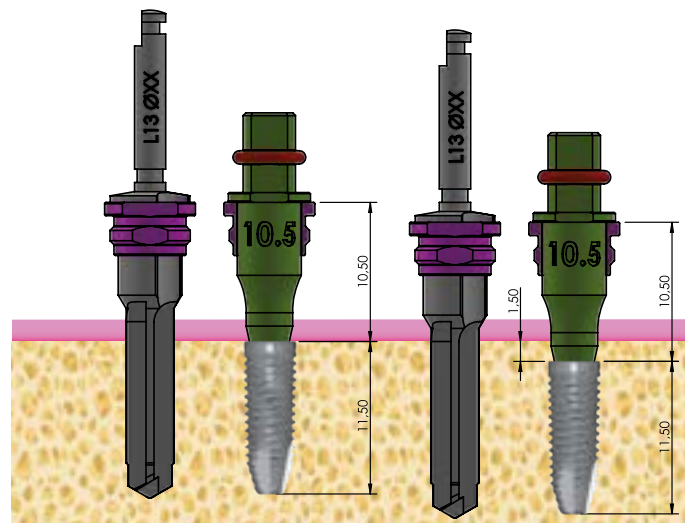


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	10.5	0 - 1.5	13

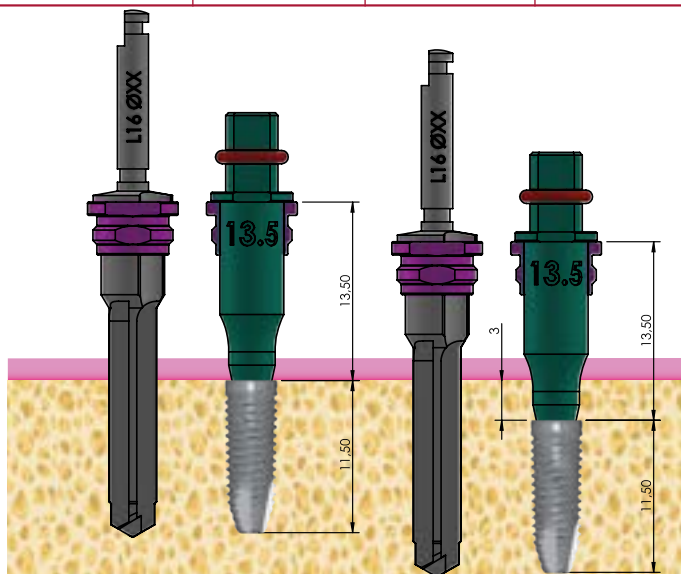


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	13.5	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 13 mm

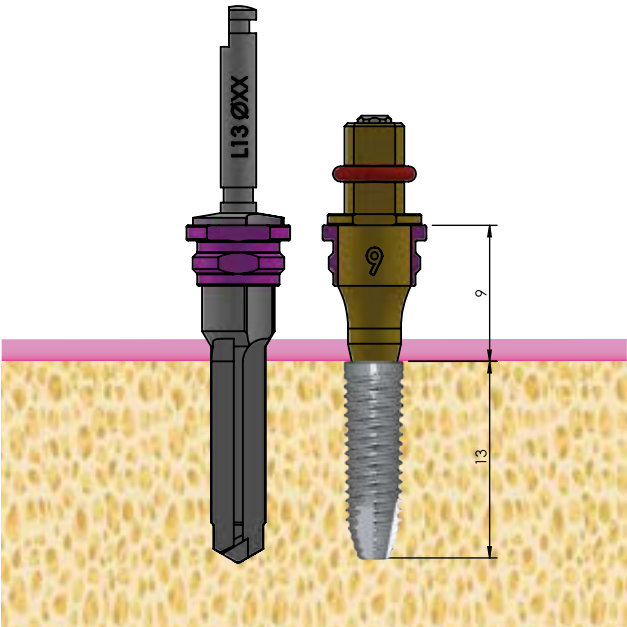


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	9	0	13

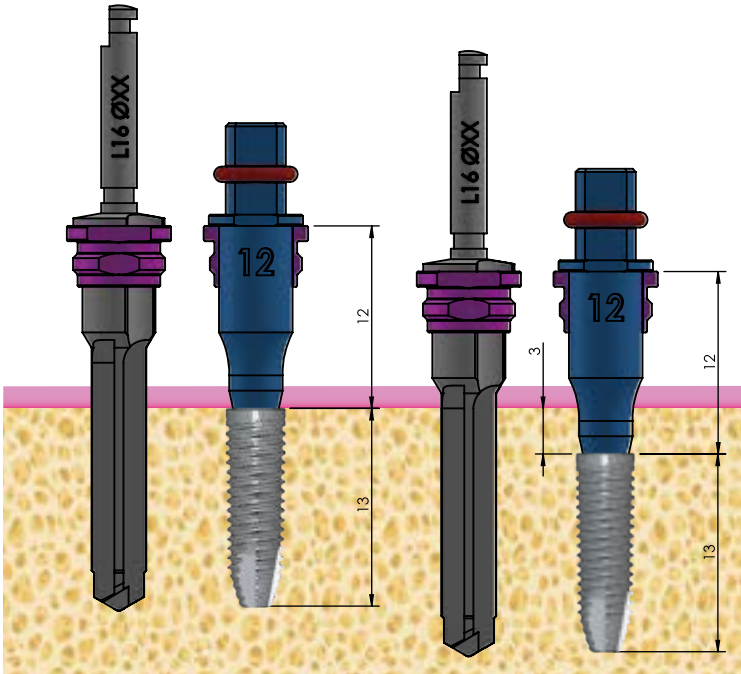

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	12	0 - 3	16

Ø XX = Ø Drill


IMPLANT LENGTH 15 mm

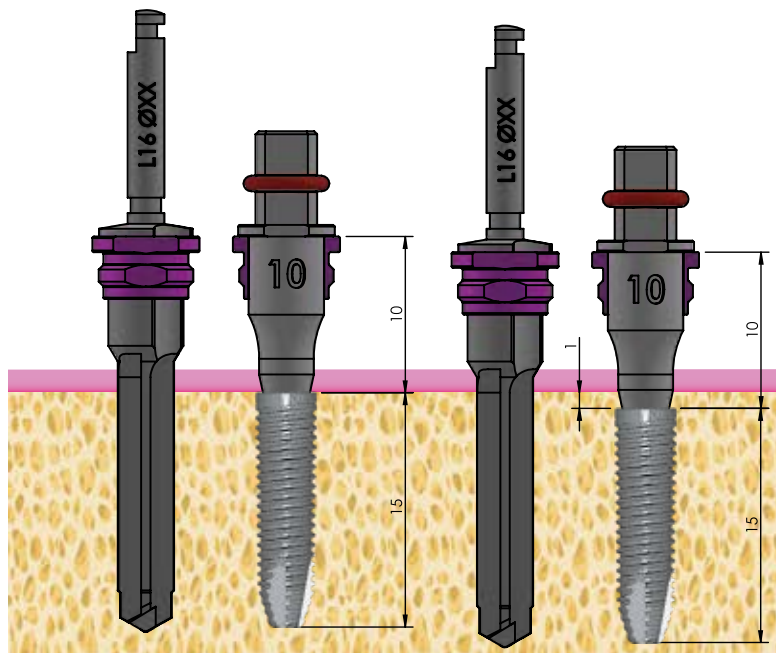


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
15	10	0 - 1	16

Ø XX = Ø Drill

IMPLANT LENGTH 16 mm

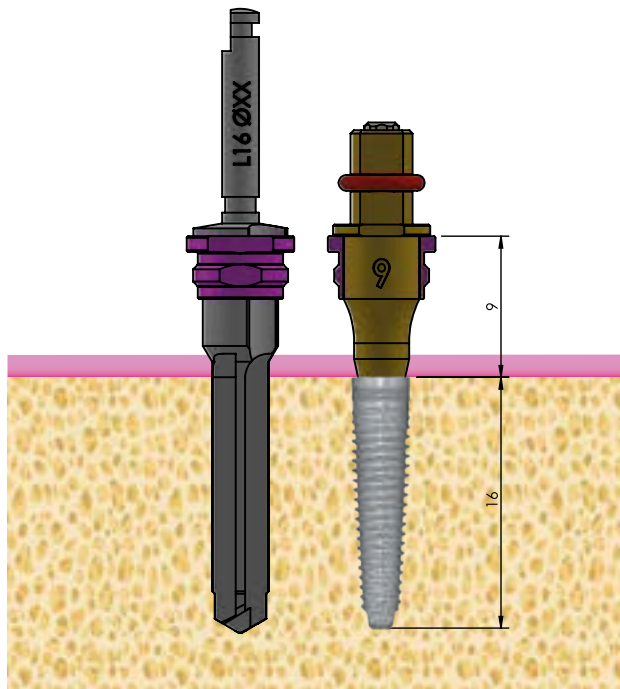


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
16	9	0	16

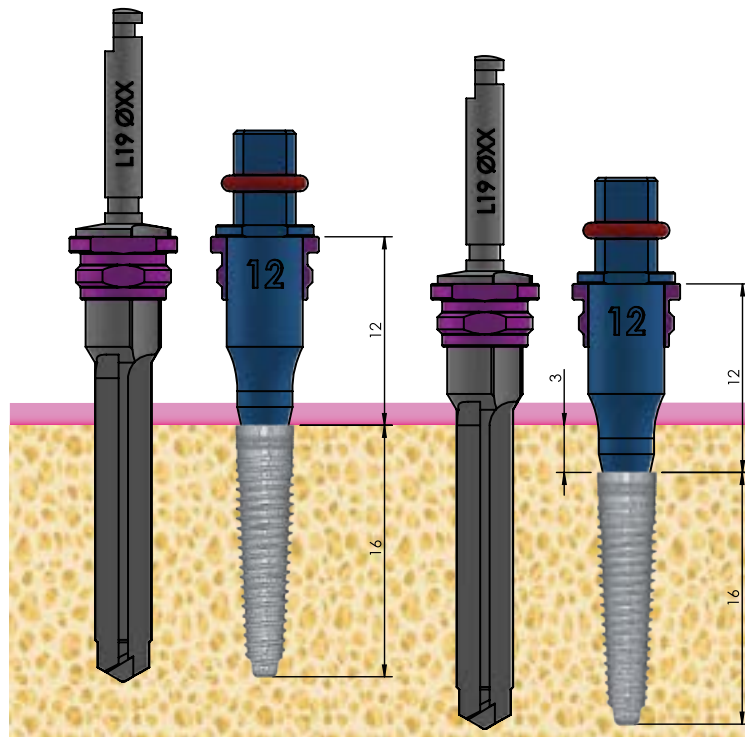

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
16	12	0 - 3	19

Ø XX = Ø Drill

IMPLANT LENGTH 18 mm

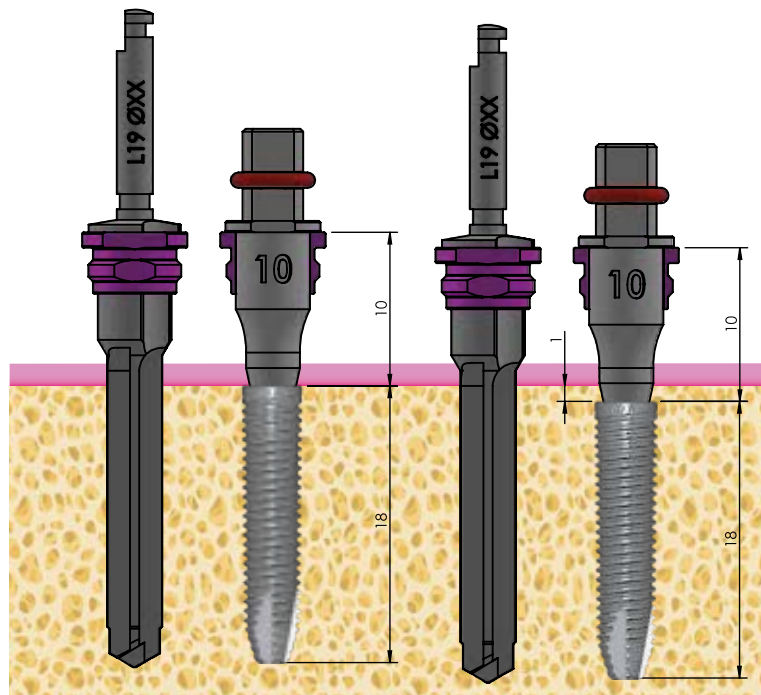

Implant Length (mm)	Abutment (mm)	Depth (mm)	Drill Length (mm)
18	10	0 - 1	19

Ø XX = Ø Drill

NOTES

[illegible]


OSSTEM - TS SYSTEM

Implant Length (mm)	Mount (mm)	Depth (mm)	Drill Length (mm)
7	10.5	0 - 1.5	8.5
	12	0 - 3	10
	13.5	0 - 3	11.5
8.5	9	0	8.5
	10.5	0 - 1.5	10
	12	0 - 3	11.5
	13.5	0 - 3	13
10	9	0	10
	10.5	0 - 1.5	11.5
	12	0 - 3	13
11.5	9	0	11.5
	10.5	0 - 1.5	13
	13.5	0 - 3	16
13	9	0	13
	12	0 - 3	16


OSSTEM TS SYSTEM	
Ø Implant (mm)	Implant Length (mm)
4/4.5	7/8.5/10/11.5/13

DAS_0030_B1


IMPLANT LENGTH 7 mm

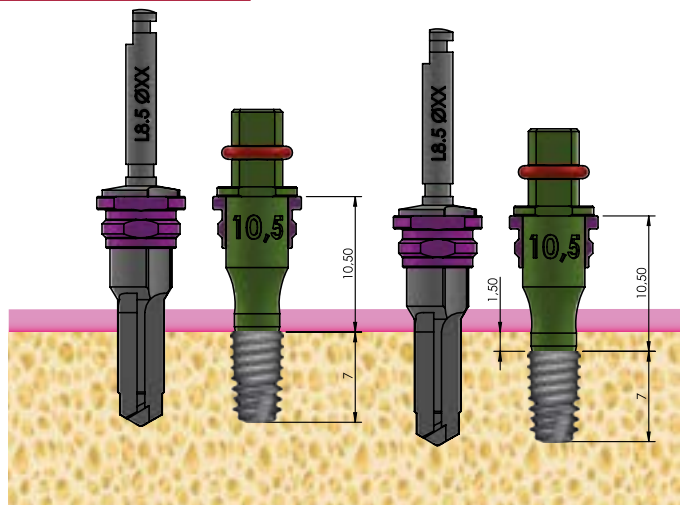

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	10.5	0 - 1.5	8.5

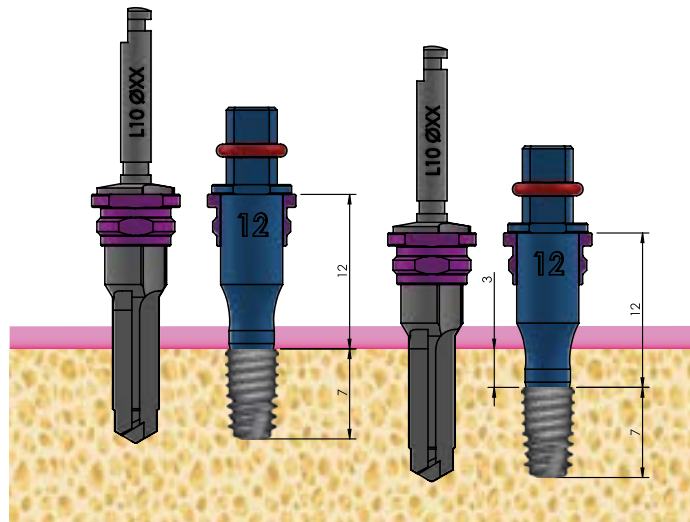


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	12	0 - 3	10

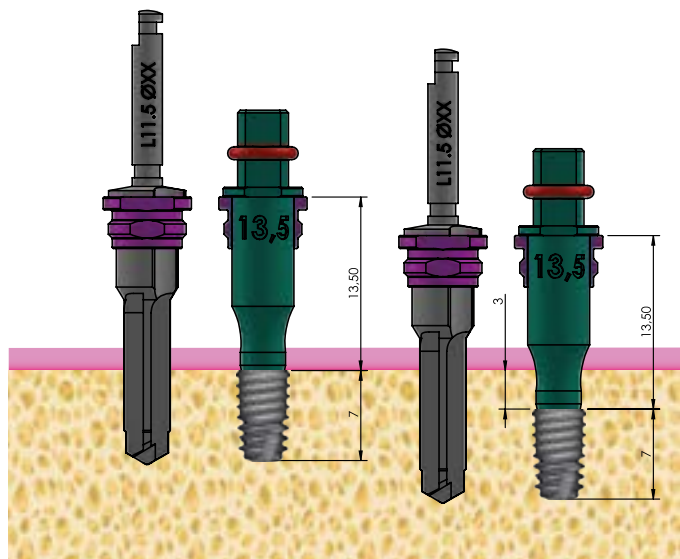


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
7	13.5	0 - 3	11.5

Ø XX = Ø Drill


IMPLANT LENGTH 8.5 mm

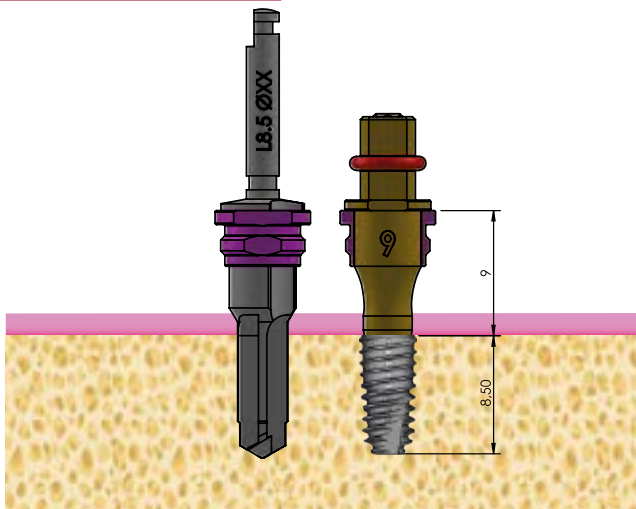


Implant Length (mm)	Moulder (mm)	Depth (mm)	Drill Length (mm)
8.5	9	0	8.5

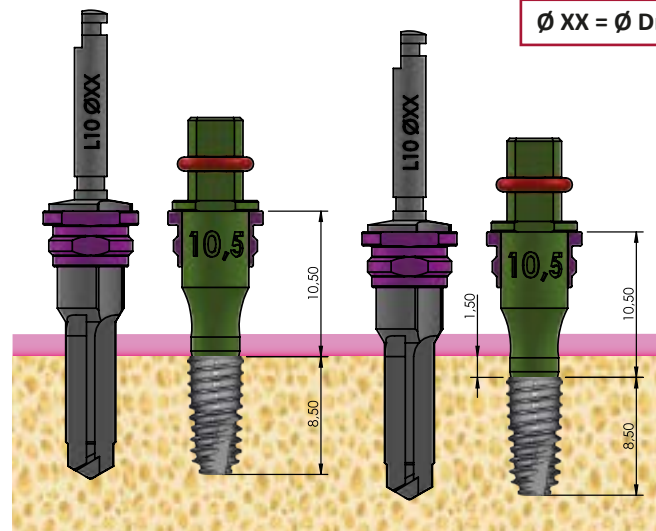 $\varnothing XX = \varnothing \text{ Drill}$ 

Implant Length (mm)	Moulder (mm)	Depth (mm)	Drill Length (mm)
8.5	10.5	0 - 1.5	10

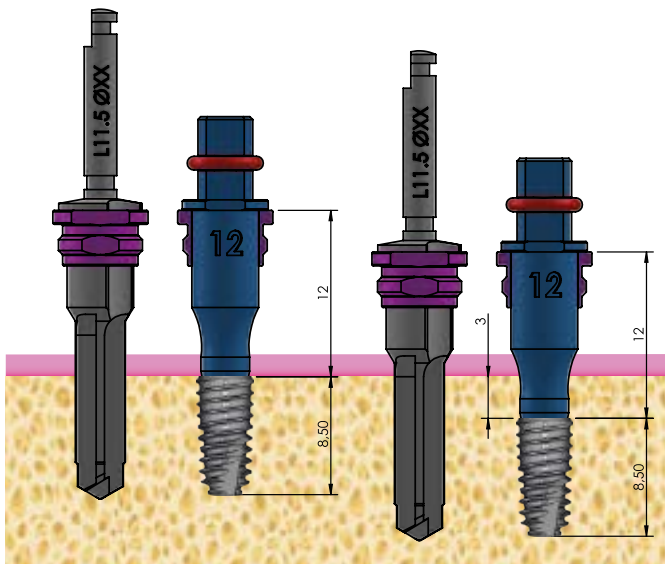


Implant Length (mm)	Moulder (mm)	Depth (mm)	Drill Length (mm)
8.5	12	0 - 3	11.5

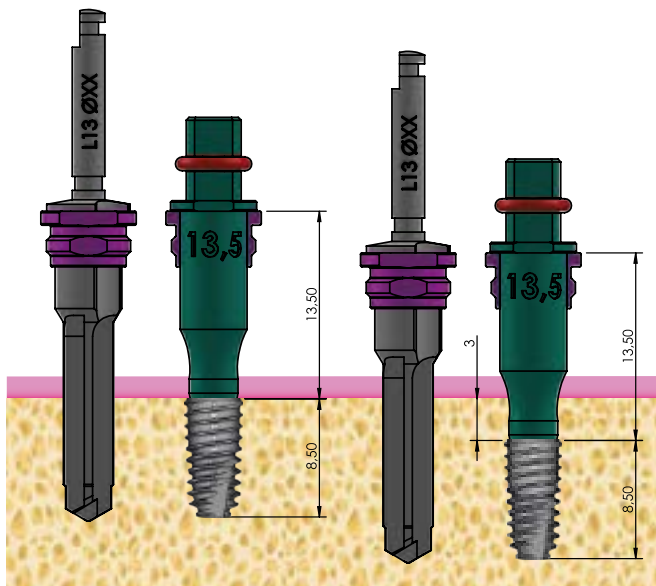


Implant Length (mm)	Moulder (mm)	Depth (mm)	Drill Length (mm)
8.5	13.5	0 - 3	13


IMPLANT LENGTH 10 mm

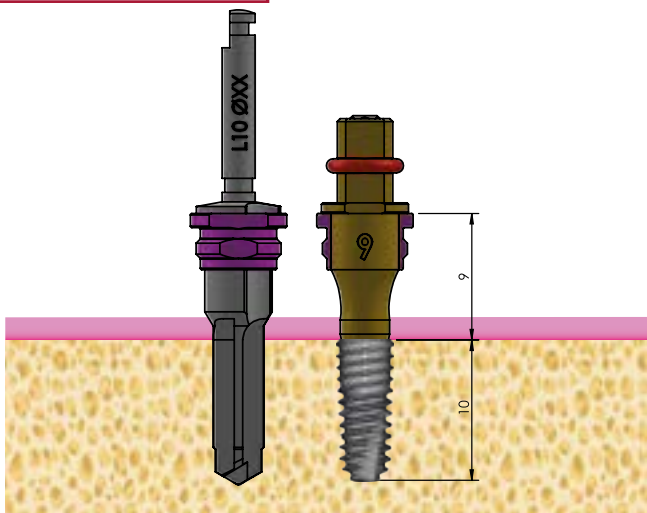

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	9	0	10

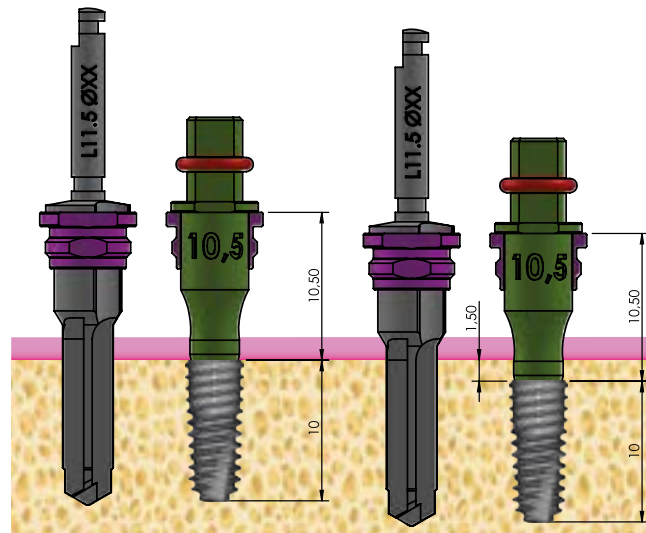


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	10.5	0 - 1.5	11.5

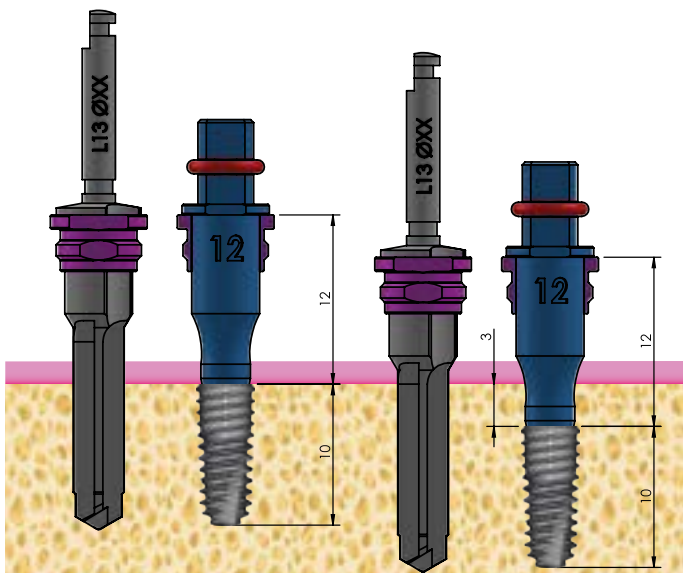


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	12	0 - 3	13

Ø XX = Ø Drill


IMPLANT LENGTH 11.5 mm

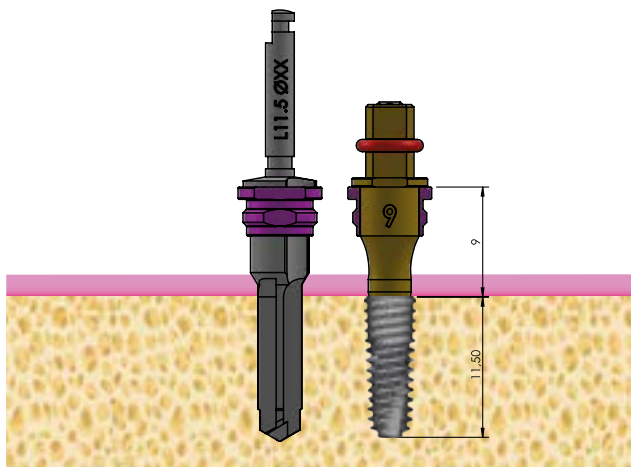

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	9	0	11.5

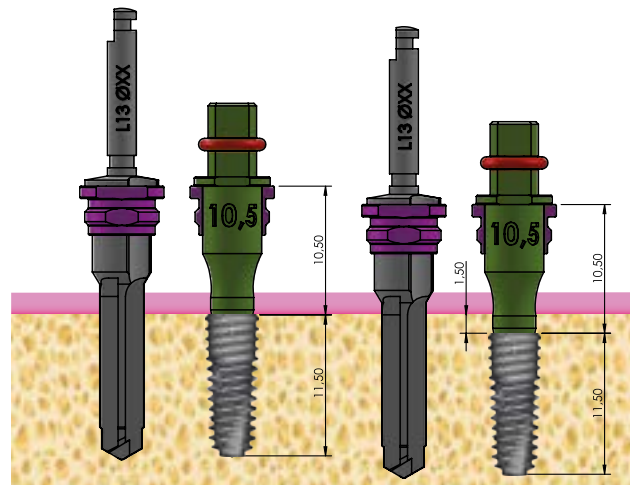


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	10.5	0 - 1.5	13

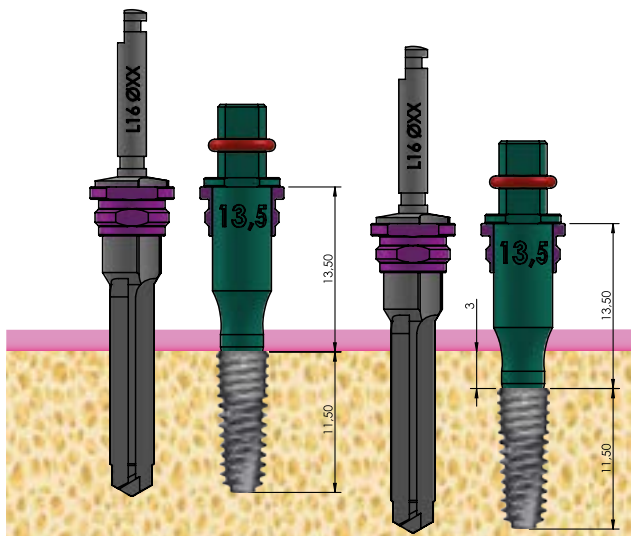


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	13.5	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 13 mm

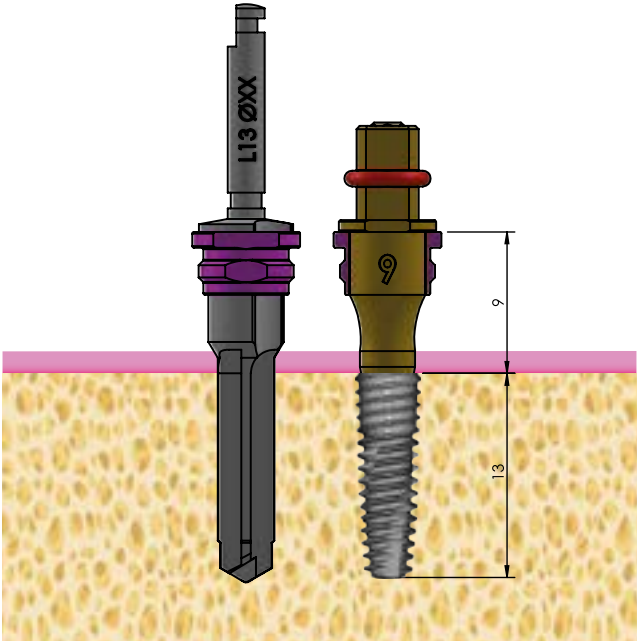

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	9	0	13

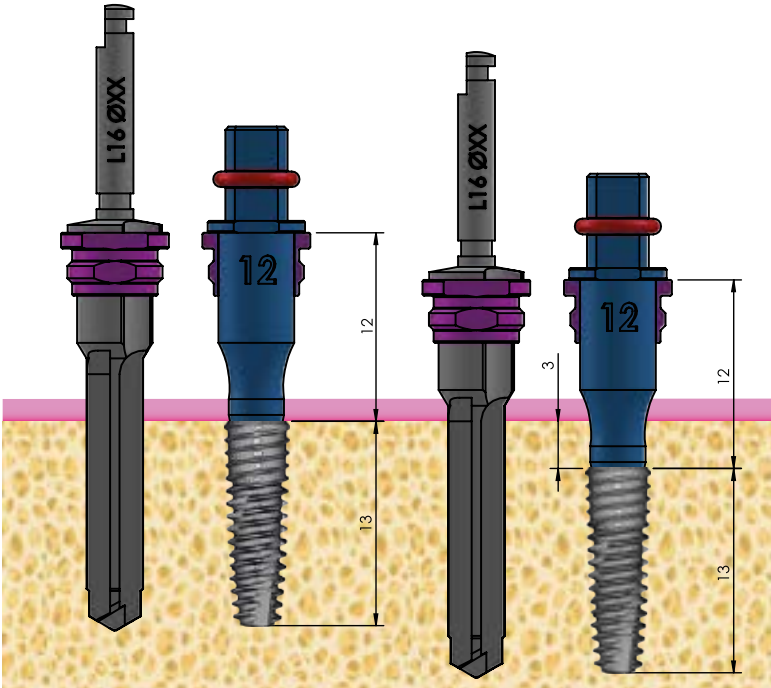

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	12	0 - 3	16

Ø XX = Ø Drill


PALTOP IH - ADVANCED/ADVANCED+/DYNAMIC/ACTIVE

Implant Length (mm)	Mount (mm)	Depth (mm)	Drill Length (mm)
6	9	0	6
8	9.5	0 - 0.5	8.5
10	9	0	10
	10.5	0 - 1.5	11.5
	12	0 - 3	13
11.5	9	0	11.5
	10.5	0 - 1.5	13
	13.5	0 - 3	16
13	9	0	13
	12	0 - 3	16
16	9	0	16
	12	0 - 3	19


PALTOP IH		
Model	Ø Implant (mm)	Implant Length (mm)
Advanced	3.75	8/10/11.5/13/16
Advanced + Dynamic	4.2	6/8/10/11.5/13/16
Active	3.75	8/10/11.5/13/16
	4.2	8/10/11.5/13/16

DAS_0040_B3

IMPLANT LENGTH 6 mm

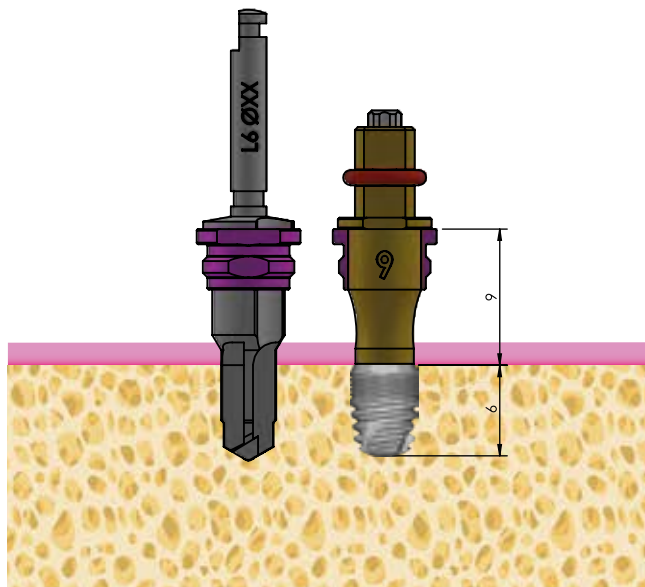

IMPLANT LENGTH 8 mm

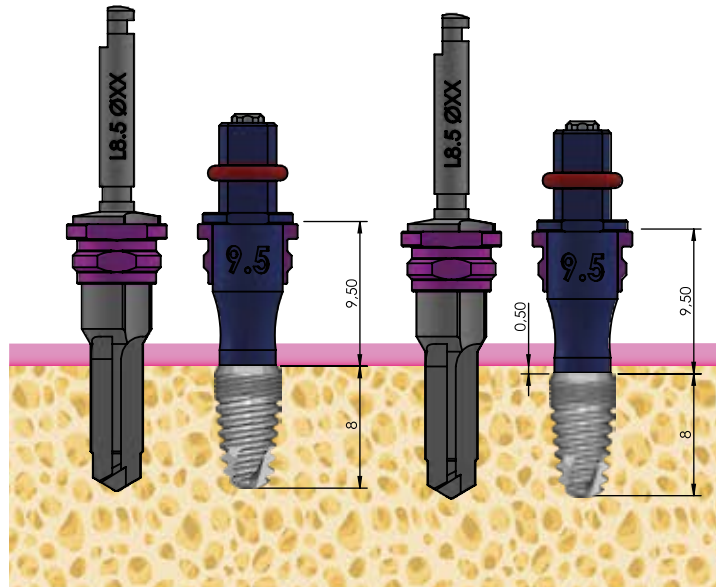


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
6	9	0	6


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
8	9.5	0 - 0.5	8.5

Ø XX = Ø Drill


IMPLANT LENGTH 10 mm

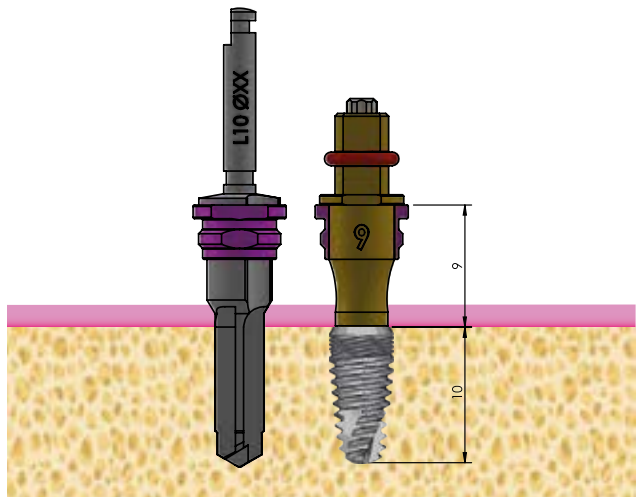

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	9	0	10

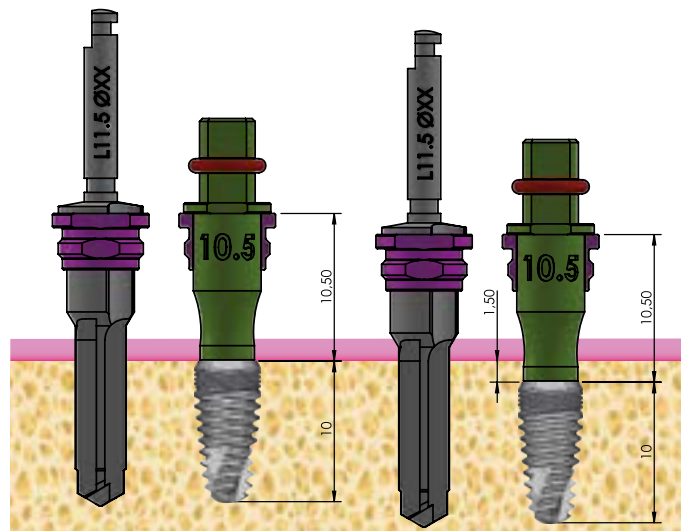


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	10.5	0 - 1.5	11.5

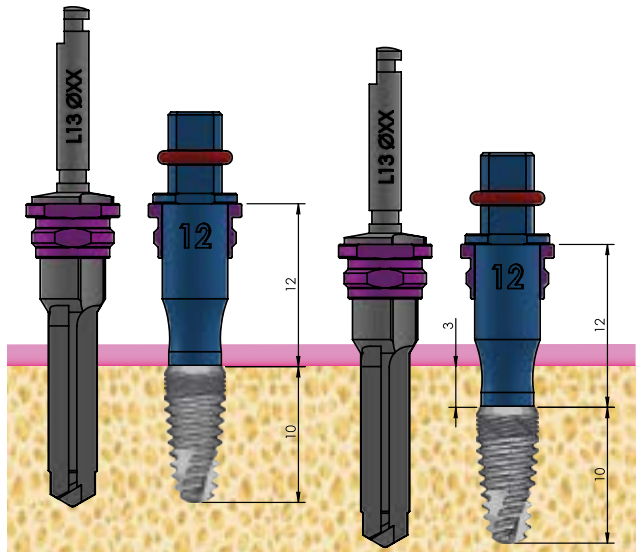


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
10	12	0 - 3	13

Ø XX = Ø Drill


IMPLANT LENGTH 11.5 mm

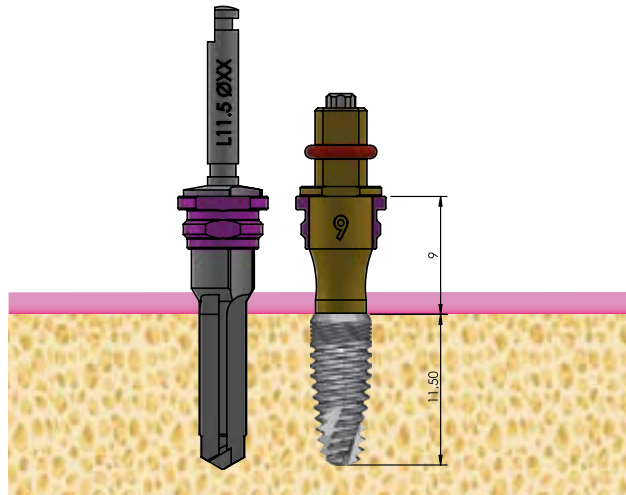

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	9	0	11.5

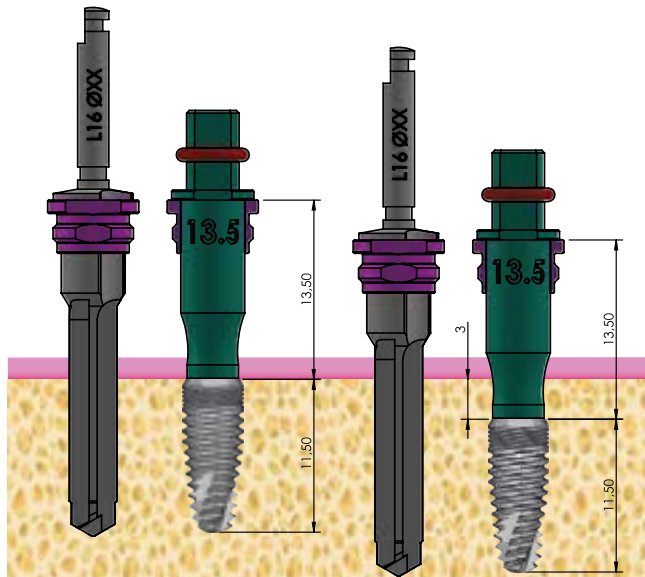


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	13.5	0 - 3	16

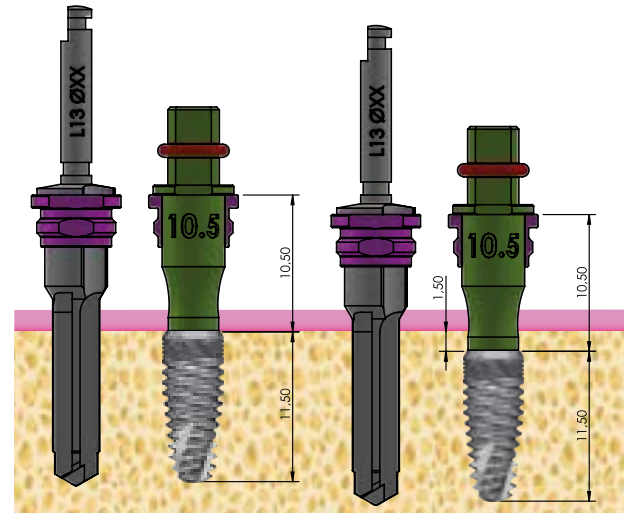


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
11.5	10.5	0 - 1.5	13

Ø XX = Ø Drill

IMPLANT LENGTH 13 mm

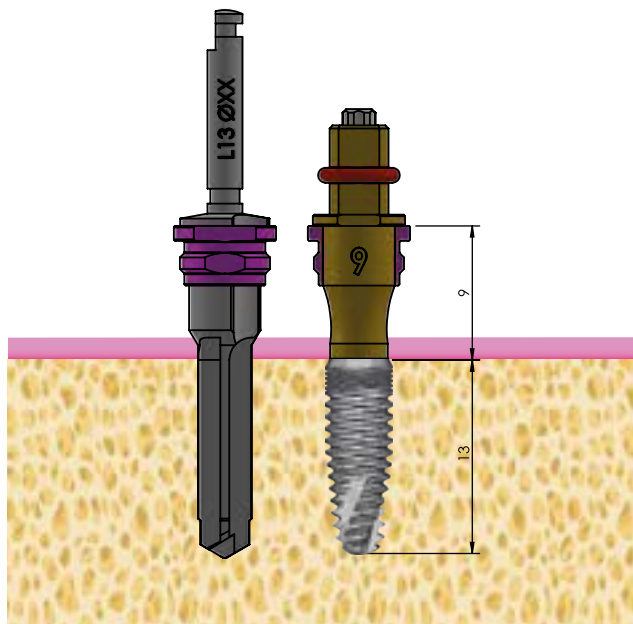


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	9	0	13

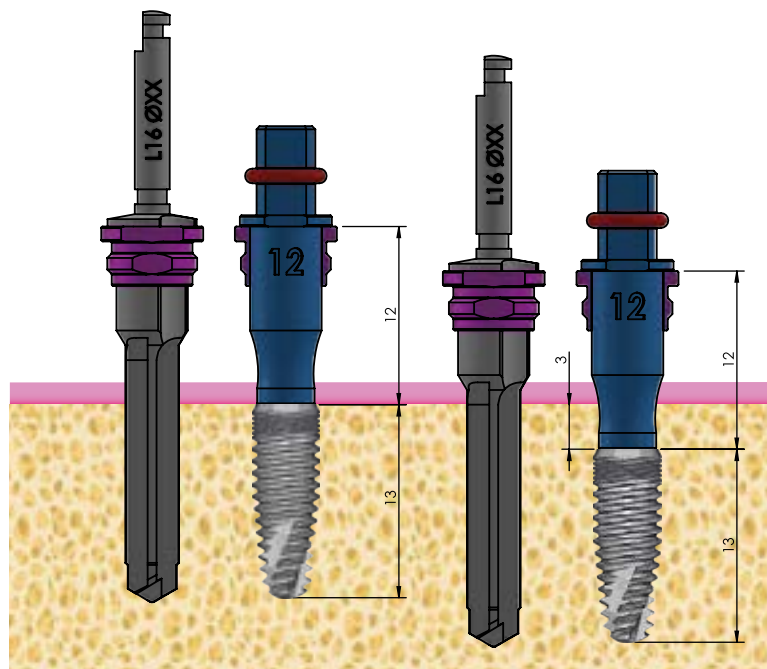


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
13	12	0 - 3	16

Ø XX = Ø Drill

IMPLANT LENGTH 16 mm

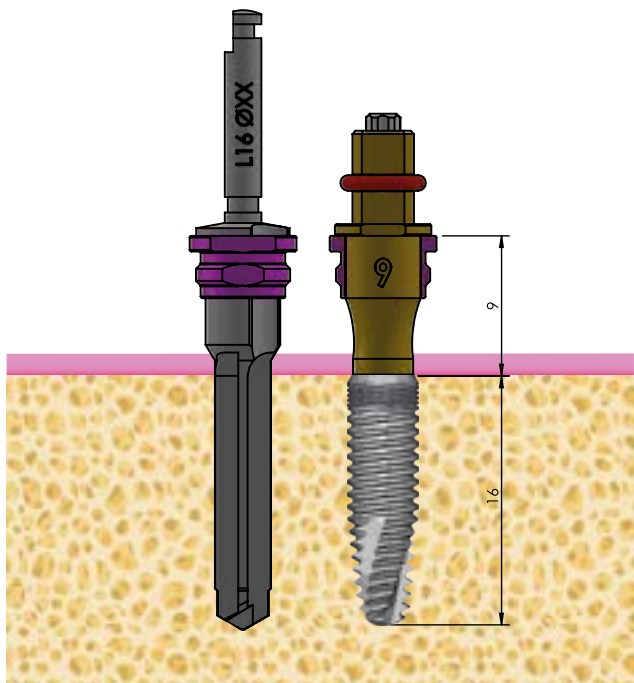


Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
16	9	0	16

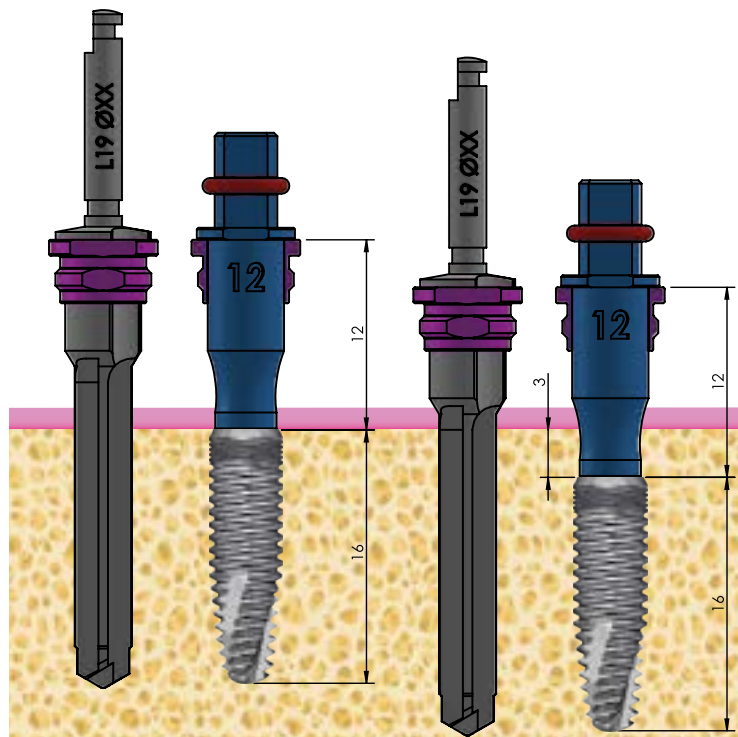

Implant Length (mm)	Mounter (mm)	Depth (mm)	Drill Length (mm)
16	12	0 - 3	19

Ø XX = Ø Drill

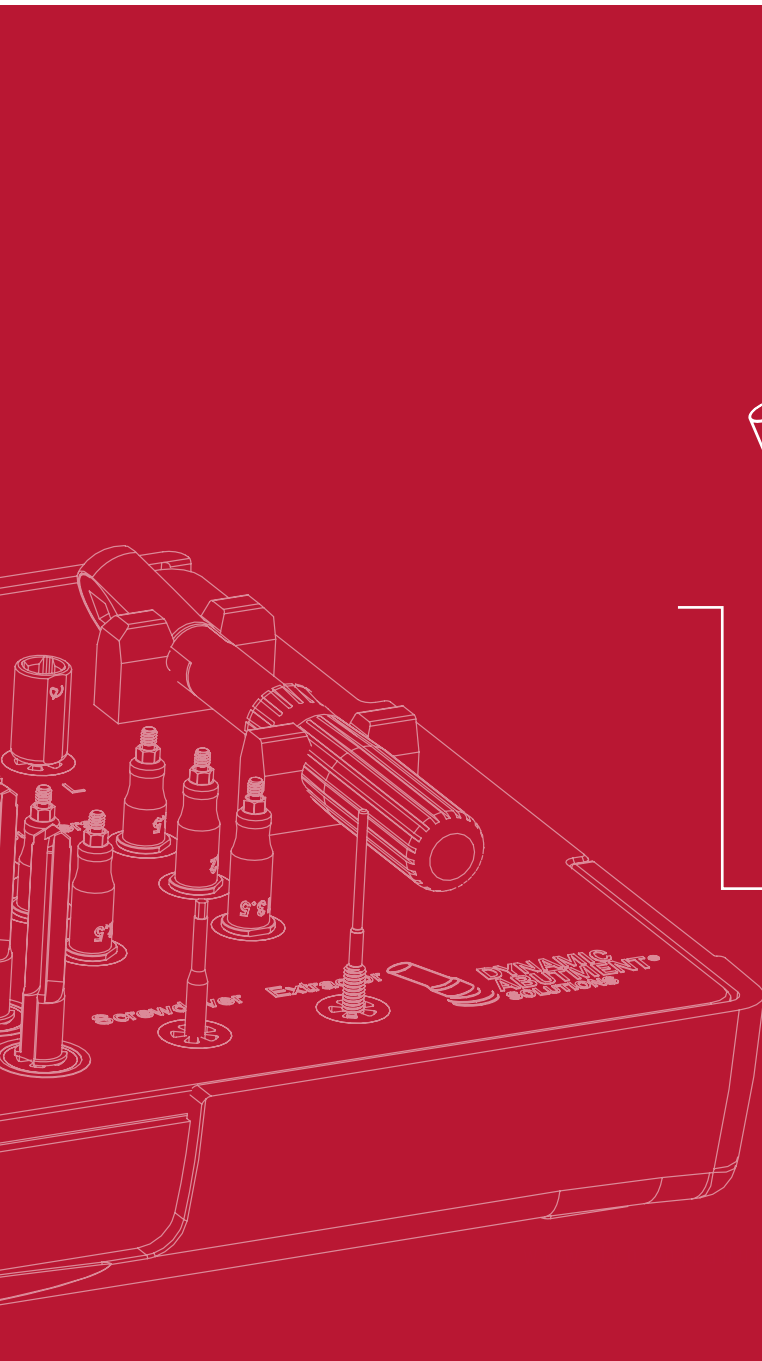

Virginia Woolf, 17
25005 · Lleida (Spain)

INTERNATIONAL +34 873 450 709
das@dynamicabutment.com

SPAIN +34 973 289 580
spain@dynamicabutment.com



[@dynamicabutmentsolutions](https://www.instagram.com/dynamicabutmentsolutions)